
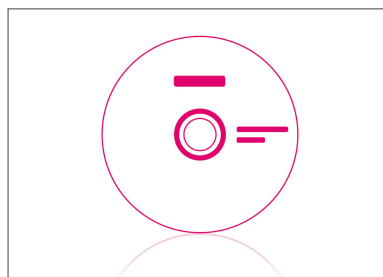
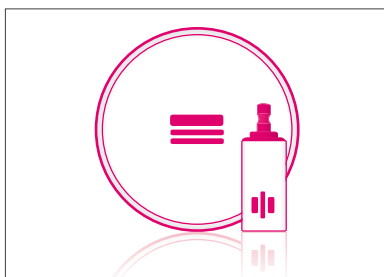
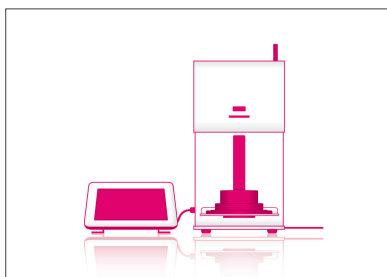
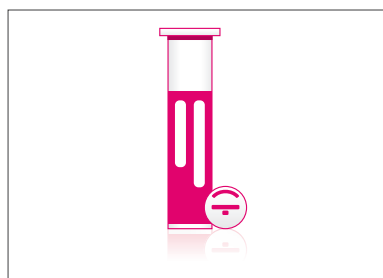
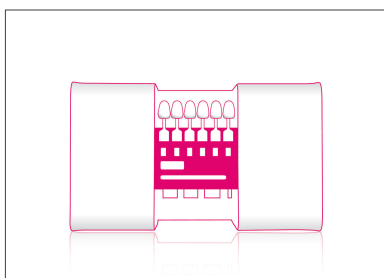


VITA Compendium September 2021

Veneering

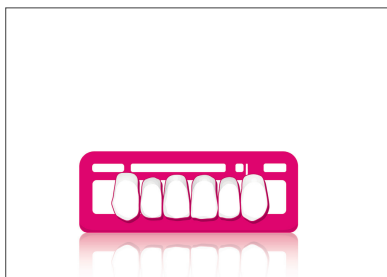

VITA – perfect match.

VITA



VITA – perfect match.

VITA

VITAVM®13

For full veneering of metal frameworks in the conventional CTE range



Product description

VITA VM 13 is a highly-esthetic, fine-structure feldspar ceramic that is perfectly adapted to the CTE value of conventional bonding alloys (approx. 13.8 – 15.2 µm/mK). Regardless of whether the conventional alloys are cast, milled or sintered.

VITA VM 13 provides a genuine alternative to all-ceramics – and is right at home wherever the masking effect of a metal framework is required.

The firing temperature (880°C, dentin firing) provides enormous processing reliability, particularly for alloys with a low solidus temperature (<1100°C).

The main components are pure-grade potash and albite feldspar materials that offer brilliant shade effects as well as unbeatable flexural strength values.

Thanks to its homogeneous, compact surface, VITA VM 13 supports outstanding grinding and polishing, particularly in situ. This ensures smooth and densely compact surface results.

Accumulation of plaque on the ceramic surface is significantly reduced, supporting the patient in caring for their high-quality restoration.

VITA VM 13 is a high-quality system component provided by the innovative VITA VM concept – the universal solution for all common indications. Metal-ceramics, all-ceramics or composite: the possibilities are endless!

VITA AKZENT Plus and VITA INTERNO stains are also available for individual finalization.

Features/Benefits

- Natural shade and light effects thanks to the fine structure
- Minimum shrinkage behavior for accurate firing
- Excellent modeling characteristics for fast and accurate application of the ceramic
- Comprehensive range of additional materials for excellent effects
- Efficient and cost-effective processing
- Can be individualized using VITA AKZENT Plus and VITA INTERNO
- Process reliability thanks to long-standing experience
- Available in VITA SYSTEM 3D-Master and VITA classical A1 – D4 shades

VITAVM®13 kits

VITAVM®13 VITA classical A1-D4®

VITAVM®13 BASIC KIT POWDER VITA classical A1-D4®

Basic kit for BASIC layering

BV13BK12CV3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE WO
16	12 g	VITA VM 13 OPAQUE A1 - D4
16	12 g	VITA VM 13 BASE DENTINE A1 - D4
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
2	12 g	VITA VM 13 ENAMEL ENL, END
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
3	12 g	VITA VM 13 CORRECTIVE COR1 - COR3
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		VITA Shade indicator for VITA VM 13 classical
1		VITA classical A1-D4 shade guide

VITAVM®13 BASIC KIT PASTE VITA classical A1-D4®

Basic kit for BASIC layering

BV13BK12PCV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
16	5 g	VITA VM 13 OPAQUE PASTE A1 - D4
16	12 g	VITA VM 13 BASE DENTINE A1 - D4
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
2	12 g	VITA VM 13 ENAMEL ENL, END
3	12 g	VITA VM 13 CORRECTIVE COR1 - COR3
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		VITA Shade indicator for VITA VM 13 classical
1		VITA classical A1-D4 shade guide

VITAVM®13 BASIC KIT SMALL POWDER VITA classical A1-D4®

Small basic kit for BASIC layering

BV13BKS12CV3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE WO
6	12 g	VITA VM 13 OPAQUE A1 - A3.5, B3, D3
6	12 g	VITA VM 13 BASE DENTINE A1 - A3.5, B3, D3
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
2	12 g	VITA VM 13 ENAMEL ENL, END
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
1	12 g	VITA VM 13 CORRECTIVE COR1
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		VITA Shade indicator for VITA VM 13 classical
1		VITA classical A1-D4 shade guide

VITAVM®13 BASIC KIT SMALL PASTE VITA classical A1-D4®

Small basic kit for BASIC layering

BV13BKS12PCV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
6	5 g	VITA VM 13 OPAQUE PASTE A1 - A3.5, B3, D3
6	12 g	VITA VM 13 BASE DENTINE A1 - A3.5, B3, D3
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
2	12 g	VITA VM 13 ENAMEL ENL, END
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
1	12 g	VITA VM 13 CORRECTIVE COR1
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		VITA Shade indicator for VITA VM 13 classical
1		VITA classical A1-D4 shade guide



VITAVM®13 BUILD UP KIT VITA classical A1-D4®

Build-up kit for BUILD UP layering

BV13BUK12CV1

Quantity	Content	Material
16	12 g	VITA VM 13 TRANSPA DENTINE A1 - D4
1	50 ml	VITA MODELLING FLUID RS

VITAVM®13 BUILD UP KIT SMALL VITA classical A1-D4®

Small build-up kit for BUILD UP layering

BV13BUKS12CV1

Quantity	Content	Material
6	12 g	VITA VM 13 TRANSPA DENTINE A1 - A3.5, B3, D3
1	50 ml	VITA MODELLING FLUID RS

VITAVM®13 CLASSICAL COLOR KIT PASTE**Add-on kit for VITA VM 13 3D-MASTER users****BV13CCK12PV1***

Quantity	Content	Material
16	5 g	VITA VM 13 OPAQUE PASTE A1 - D4
16	12 g	VITA VM 13 BASE DENTINE A1 - D4
16	12 g	VITA VM 13 TRANSPA DENTINE A1 - D4
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1		VITA Shade indicator for VITA VM 13 classical
1		VITA classical A1-D4 shade guide

* Assortment for VITA VM 13 SYSTEM 3D-MASTER customers who would like to add VITA classical A1-D4 to their existing range of materials

**VITAVM®13 CLASSICAL COLOR KIT POWDER****Add-on kit for VITA VM 13 3D-MASTER users****BV13CCK12V1***

Quantity	Content	Material
16	12 g	VITA VM 13 OPAQUE A1 - D4
16	12 g	VITA VM 13 BASE DENTINE A1 - D4
16	12 g	VITA VM 13 TRANSPA DENTINE A1 - D4
1	50 ml	VITA VM OPAQUE FLUID
1	50 ml	VITA MODELLING FLUID RS
1		VITA Shade indicator for VITA VM 13 classical
1		VITA classical A1-D4 shade guide

* Assortment for VITA VM 13 SYSTEM 3D-MASTER customers who would like to add VITA classical A1-D4 to their existing range of materials

VITAVM®13 ONE COLOR KIT A2 PASTE

One-shade kit for initial tests

BV130CKA2PV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
1	5 g	VITA VM 13 OPAQUE PASTE A2
1	12 g	VITA VM 13 CHROMA PLUS CP2
1	12 g	VITA VM 13 BASE DENTINE A2
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	4 g	VITA AKZENT Plus GLAZE PASTE
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1	15 ml	VITA AKZENT Plus PASTE FLUID
1		VITA shade indicator stick A2



VITAVM®13 ONE COLOR KIT A2 POWDER

One-shade kit for initial tests

BV130CKA2V3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE W0
1	12 g	VITA VM 13 OPAQUE A2
1	12 g	VITA VM 13 CHROMA PLUS CP2
1	12 g	VITA VM 13 BASE DENTINE A2
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1	20 ml	VITA AKZENT Plus POWDER FLUID
1		VITA shade indicator stick A2

VITAVM®13 ONE COLOR KIT A3 PASTE

One-shade kit for initial tests

BV130CKA3PV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
1	5 g	VITA VM 13 OPAQUE PASTE A3
2	12 g	VITA VM 13 CHROMA PLUS CP2 - CP3
1	12 g	VITA VM 13 BASE DENTINE A3
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	4 g	VITA AKZENT Plus GLAZE PASTE
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1	15 ml	VITA AKZENT Plus PASTE FLUID
1		VITA shade indicator stick A3

VITAVM®13 ONE COLOR KIT A3 POWDER

One-shade kit for initial tests

BV130CKA3V3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE W0
1	12 g	VITA VM 13 OPAQUE A3
2	12 g	VITA VM 13 CHROMA PLUS CP2 - CP3
1	12 g	VITA VM 13 BASE DENTINE A3
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1	20 ml	VITA AKZENT Plus POWDER FLUID
1		VITA shade indicator stick A3

VITAVM®13 VITA SYSTEM 3D-MASTER®
VITAVM®13 BASIC KIT PASTE VITA SYSTEM 3D-MASTER®
Basic kit for BASIC layering
BV13BK12PV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
5	5 g	VITA VM 13 OPAQUE PASTE OP1 - OP5
3	5 g	VITA VM 13 SUN OPAQUE PASTE SO1 - SO3
26	12 g	VITA VM 13 BASE DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
2	12 g	VITA VM 13 ENAMEL ENL, END
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
3	12 g	VITA VM 13 CORRECTIVE COR1 - COR3
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 13 3D-MASTER
1		VITA Toothguide 3D-Master with Bleached Shade Guide

VITAVM®13 BASIC KIT POWDER VITA SYSTEM 3D-MASTER®

Basic kit for BASIC layering

BV13BK12V3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE WO
5	12 g	VITA VM 13 OPAQUE OP1 - OP5
3	12 g	VITA VM 13 SUN OPAQUE SO1 - SO3
26	12 g	VITA VM 13 BASE DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
2	12 g	VITA VM 13 ENAMEL ENL, END
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
3	12 g	VITA VM 13 CORRECTIVE COR1 - COR3
1		VITA Modisol®
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 13 3D-MASTER
1		VITA Toothguide 3D-Master with Bleached Shade Guide

VITAVM®13 BASIC KIT SMALL PASTE VITA SYSTEM 3D-MASTER®

Small basic kit for BASIC layering

BV13BKS12PV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
4	5 g	VITA VM 13 OPAQUE PASTE OP1 - OP4
3	5 g	VITA VM 13 SUN OPAQUE PASTE SO1 - SO3
15	12 g	VITA VM 13 BASE DENTINE 1M1 - 4M3
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
2	12 g	VITA VM 13 ENAMEL ENL, END
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
1	12 g	VITA VM 13 CORRECTIVE COR1
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 13 3D-MASTER
1		VITA Toothguide 3D-Master with Bleached Shade Guide

VITAVM®13 BASIC KIT SMALL POWDER VITA SYSTEM 3D-MASTER®
Small basic kit for BASIC layering
BV13BKS12V3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE WO
4	12 g	VITA VM 13 OPAQUE OP1 - OP4
3	12 g	VITA VM 13 SUN OPAQUE SO1 - SO3
15	12 g	VITA VM 13 BASE DENTINE 1M1 - 4M3
3	12 g	VITA VM 13 SUN DENTINE SD1 - SD3
2	12 g	VITA VM 13 ENAMEL ENL, END
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
1	12 g	VITA VM 13 CORRECTIVE COR1
5	12 g	VITA VM 13 CHROMA PLUS CP1 - CP5
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 13 3D-MASTER
1		VITA Toothguide 3D-Master with Bleached Shade Guide

VITAVM®13 BLEACHED COLOR KIT PASTE VITA SYSTEM 3D-MASTER®
Very light shades to reproduce bleached teeth
BV13BLKPV1

Quantity	Content	Material
1	5 g	VITA VM 13 OPAQUE PASTE OP0
3	12 g	VITA VM 13 BASE DENTINE 0M1 - 0M3
3	12 g	VITA VM 13 TRANSPA DENTINE 0M1 - 0M3
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
1	15 ml	VITA VM PASTE FLUID
1	50 ml	VITA MODELLING FLUID RS
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA BLEACHED SHADE GUIDE, Shade Group 0



VITAVM®13 BLEACHED COLOR KIT POWDER VITA SYSTEM 3D-MASTER®
Very light shades to reproduce bleached teeth
BV13BLKV1

Quantity	Content	Material
1	12 g	VITA VM 13 OPAQUE OP0
3	12 g	VITA VM 13 BASE DENTINE 0M1 - 0M3
3	12 g	VITA VM 13 TRANSPA DENTINE 0M1 - 0M3
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 NEUTRAL NT
1	12 g	VITA VM 13 WINDOW WIN
1	50 ml	VITA VM OPAQUE FLUID
1	50 ml	VITA MODELLING FLUID RS
1		VITA BLEACHED SHADE GUIDE, Shade Group 0



VITAVM®13 BUILD UP KIT VITA SYSTEM 3D-MASTER®
Build-up kit for BUILD UP layering
BV13BUK12V1

Quantity	Content	Material
26	12 g	VITA VM 13 TRANSPA DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
1	50 ml	VITA MODELLING FLUID RS

VITAVM®13 BUILD UP KIT SMALL VITA SYSTEM 3D-MASTER®
Small build-up kit for BUILD UP layering
BV13BUKS12V1

Quantity	Content	Material
15	12 g	VITA VM 13 TRANSPA DENTINE 1M1 - 4M3
1	50 ml	VITA MODELLING FLUID RS

VITAVM®13 ONE COLOR KIT 2M2 PASTE

One-shade kit for initial tests

BV130CK2M2PV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
1	5 g	VITA VM 13 OPAQUE PASTE OP2
1	12 g	VITA VM 13 BASE DENTINE 2M2
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	4 g	VITA AKZENT Plus GLAZE PASTE
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1	15 ml	VITA AKZENT Plus PASTE FLUID
1		VITA shade indicator stick 2M2



VITAVM®13 ONE COLOR KIT 2M2 POWDER

One-shade kit for initial tests

BV130CK2M2V3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE W0
1	12 g	VITA VM 13 OPAQUE OP2
1	12 g	VITA VM 13 BASE DENTINE 2M2
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1	20 ml	VITA AKZENT Plus POWDER FLUID
1		VITA shade indicator stick 2M2

VITAVM®13 ONE COLOR KIT 3M2 PASTE

One-shade kit for initial tests

BV130CK3M2PV3

Quantity	Content	Material
1	7 g	VITA VM 13 WASH OPAQUE PASTE W0
1	5 g	VITA VM 13 OPAQUE PASTE OP3
1	12 g	VITA VM 13 BASE DENTINE 3M2
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	4 g	VITA AKZENT Plus GLAZE PASTE
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA VM PASTE FLUID
1	15 ml	VITA AKZENT Plus PASTE FLUID
1		VITA shade indicator stick 3M2

VITAVM®13 ONE COLOR KIT 3M2 POWDER

One-shade kit for initial tests

BV130CK3M2V3

Quantity	Content	Material
1	12 g	VITA VM 13 WASH OPAQUE W0
1	12 g	VITA VM 13 OPAQUE OP3
1	12 g	VITA VM 13 BASE DENTINE 3M2
1	12 g	VITA VM 13 ENAMEL ENL
1	12 g	VITA VM 13 WINDOW WIN
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA VM OPAQUE FLUID
1	20 ml	VITA AKZENT Plus POWDER FLUID
1		VITA shade indicator stick 3M2

VITAVM®13 additional powder

compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades

VITAVM®13 GINGIVA KIT PASTE

Natural-looking gingiva pastes

BV13GKPV1

Quantity	Content	Material
5	12 g	VITA VM 13 GINGIVA G1 - G5
2	5 g	VITA VM 13 GINGIVA OPAQUE PASTE GOL, GOD
1		Shade guide for VITA VM GINGIVA / MARGIN



VITAVM®13 GINGIVA KIT
Natural-looking gingiva pastes
BV13GKV1

Quantity	Content	Material
5	12 g	VITA VM 13 GINGIVA G1 - G5
2	12 g	VITA VM 13 GINGIVA OPAQUE GOL, GOD
1		Shade guide for VITA VM GINGIVA / MARGIN



VITAVM®13 MARGIN KIT
Ceramic for creating the transition between crown and prepared tooth
BV13MKV1

Quantity	Content	Material
6	12 g	VITA VM 13 MARGIN M1, M4 - M5, M7 - M9
1		Shade guide for VITA VM GINGIVA / MARGIN



VITAVM®13 PROFESSIONAL KIT SMALL
Additional materials for achieving special shade effects and individualizing the restoration
BV13PKS

Quantity	Content	Material
5	12 g	VITA VM 13 EFFECT CHROMA EC1, EC4, EC6, EC8 - EC9
1	12 g	VITA VM 13 MAMELON MM2
1	12 g	VITA VM 13 EFFECT PEARL EP1
1	12 g	VITA VM 13 EFFECT OPAL EO2
7	12 g	VITA VM 13 EFFECT ENAMEL EE1, EE3, EE7 - EE11
1		Shade guide for VITA VM PROFESSIONAL KIT SMALL



VITAVM®13 PROFESSIONAL KIT
Additional materials for achieving special shade effects and individualizing the restoration
BV13PKV2

Quantity	Content	Material
11	12 g	VITA VM 13 EFFECT CHROMA EC1 - EC11
3	12 g	VITA VM 13 MAMELON MM1 - MM3
3	12 g	VITA VM 13 EFFECT PEARL EP1 - EP3
5	12 g	VITA VM 13 EFFECT OPAL EO1 - EO5
11	12 g	VITA VM 13 EFFECT ENAMEL EE1 - EE11
6	12 g	VITA VM 13 EFFECT LINER EL1 - EL6
1		Shade guide for VITA VM EFFECT CHROMA
1		Shade guide for VITA VM EFFECT LINER / EFFECT OPAL / MAMELON
1		Shade guide for VITA VM EFFECT ENAMEL / EFFECT PEARL

VITAVM®13 products**VITAVM®13 VITA classical A1-D4®**

VITA VM 13 OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
A1	B453015
A2	B453025
A3	B453035
A3.5	B453045
A4	B453055
B1	B453065
B2	B453075
B3	B453085
B4	B453095
C1	B453105
C2	B453115
C3	B453125
C4	B453135
D2	B453145
D3	B453155
D4	B453165

VITA VM 13 OPAQUE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
A1	B4540112	B4540150
A2	B4540212	B4540250
A3	B4540312	B4540350
A3.5	B4540412	B4540450
A4	B4540512	B4540550
B1	B4540612	B4540650
B2	B4540712	B4540750
B3	B4540812	B4540850
B4	B4540912	B4540950
C1	B4541012	B4541050
C2	B4541112	B4541150
C3	B4541212	B4541250
C4	B4541312	B4541350
D2	B4541412	B4541450
D3	B4541512	B4541550
D4	B4541612	B4541650

VITA VM 13 BASE DENTINE, 12 g / 50 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
A1	ENL	B4543112	B4543150
A2	ENL	B4543212	B4543250
A3	ENL	B4543312	B4543350
A3.5	END	B4543412	B4543450
A4	END	B4543512	B4543550
B1	END	B4543612	B4543650
B2	END	B4543712	B4543750
B3	END	B4543812	B4543850
B4	END	B4543912	B4543950
C1	END	B4544012	B4544050
C2	END	B4544112	B4544150
C3	ENL	B4544212	B4544250
C4	ENL	B4544312	B4544350
D2	END	B4544412	B4544450
D3	END	B4544512	B4544550
D4	END	B4544612	B4544650

VITA VM 13 TRANSPA DENTINE, 12 g / 50 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
A1	ENL	B4545112	B4545150
A2	ENL	B4545212	B4545250
A3	ENL	B4545312	B4545350
A3.5	END	B4545412	B4545450
A4	END	B4545512	B4545550
B1	END	B4545612	B4545650
B2	END	B4545712	B4545750
B3	END	B4545812	B4545850
B4	END	B4545912	B4545950
C1	END	B4546012	B4546050
C2	END	B4546112	B4546150
C3	ENL	B4546212	B4546250
C4	ENL	B4546312	B4546350
D2	END	B4546412	B4546450
D3	END	B4546512	B4546550
D4	END	B4546612	B4546650

VITAVM® 13 VITA SYSTEM 3D-MASTER®

VITA VM 13 BASE DENTINE, 12 g / 50 g / 250 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
0M1	ENL	B4503112	B4503150	-
0M2	ENL	B4503212	B4503250	-
0M3	ENL	B4503312	B4503350	-
1M1	ENL	B4503412	B4503450	-
1M2	ENL	B4503512	B4503550	-
2L1.5	ENL	B4503612	B4503650	-
2L2.5	ENL	B4503712	B4503750	-
2M1	ENL	B4503812	B4503850	-
2M2	ENL	B4503912	B4503950	B45039250
2M3	ENL	B4504012	B4504050	B45040250
2R1.5	ENL	B4504112	B4504150	-
2R2.5	ENL	B4504212	B4504250	-
3L1.5	ENL	B4504312	B4504350	-
3L2.5	ENL	B4504412	B4504450	-
3M1	ENL	B4504512	B4504550	-
3M2	ENL	B4504612	B4504650	B45046250
3M3	ENL	B4504712	B4504750	B45047250

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
3R1.5	ENL	B4504812	B4504850	-
3R2.5	ENL	B4504912	B4504950	-
4L1.5	END	B4505012	B4505050	-
4L2.5	END	B4505112	B4505150	-
4M1	END	B4505212	B4505250	-
4M2	END	B4505312	B4505350	-
4M3	END	B4505412	B4505450	-
4R1.5	END	B4505512	B4505550	-
4R2.5	END	B4505612	B4505650	-
5M1	END	B4505712	B4505750	-
5M2	END	B4505812	B4505850	-
5M3	END	B4505912	B4505950	-

VITA VM 13 TRANSPA DENTINE, 12 g / 50 g / 250 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
0M1	ENL	B4506112	B4506150	-
0M2	ENL	B4506212	B4506250	-
0M3	ENL	B4506312	B4506350	-
1M1	ENL	B4506412	B4506450	-
1M2	ENL	B4506512	B4506550	-
2L1.5	ENL	B4506612	B4506650	-
2L2.5	ENL	B4506712	B4506750	-
2M1	ENL	B4506812	B4506850	-
2M2	ENL	B4506912	B4506950	B45069250
2M3	ENL	B4507012	B4507050	B45070250
2R1.5	ENL	B4507112	B4507150	-
2R2.5	ENL	B4507212	B4507250	-
3L1.5	ENL	B4507312	B4507350	-
3L2.5	ENL	B4507412	B4507450	-
3M1	ENL	B4507512	B4507550	-
3M2	ENL	B4507612	B4507650	B45076250
3M3	ENL	B4507712	B4507750	B45077250
3R1.5	ENL	B4507812	B4507850	-
3R2.5	ENL	B4507912	B4507950	-
4L1.5	END	B4508012	B4508050	-
4L2.5	END	B4508112	B4508150	-
4M1	END	B4508212	B4508250	-
4M2	END	B4508312	B4508350	-
4M3	END	B4508412	B4508450	-

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
4R1.5	END	B4508512	B4508550	-
4R2.5	END	B4508612	B4508650	-
5M1	END	B4508712	B4508750	-
5M2	END	B4508812	B4508850	-
5M3	END	B4508912	B4508950	-

VITA VM 13 OPAQUE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
OP0	B4525012	-
OP1	B4525112	B4525150
OP2	B4525212	B4525250
OP3	B4525312	B4525350
OP4	B4525412	B4525450
OP5	B4525512	-

VITA VM 13 OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
OP0	B453505
OP1	B453515
OP2	B453525
OP3	B453535
OP4	B453545
OP5	B453555

VITAVM®13 additional powder
compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades

VITA VM 13 EFFECT ENAMEL, 12 g

Description	Prod. no. 12 g
EE1	B4509112
EE2	B4509212
EE3	B4509312
EE4	B4509412
EE5	B4509512
EE6	B4509612
EE7	B4509712
EE8	B4509812
EE9	B4509912
EE10	B4510012
EE11	B4510112

VITA VM 13 EFFECT CHROMA, 12 g

Description	Prod. no. 12 g
EC1	B4512112
EC2	B4512212
EC3	B4512312
EC4	B4512412
EC5	B4512512
EC6	B4512612
EC7	B4512712
EC8	B4512812
EC9	B4512912
EC10	B4513012
EC11	B4513112

VITA VM 13 EFFECT LINER, 12 g

Description	Prod. no. 12 g
EL1	B4514112
EL2	B4514212
EL3	B4514312
EL4	B4514412
EL5	B4514512
EL6	B4514612

VITA VM 13 EFFECT PEARL, 12 g

Description	Prod. no. 12 g
EP1	B4515112
EP2	B4515212
EP3	B4515312

VITA VM 13 MAMELON, 12 g

Description	Prod. no. 12 g
MM1	B4516112
MM2	B4516212
MM3	B4516312

VITA VM 13 EFFECT OPAL, 12 g

Description	Prod. no. 12 g
E01	B4517112
E02	B4517212
E03	B4517312
E04	B4517412
E05	B4517512

VITA VM 13 WINDOW, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
WIN	B4518112	B4518150	B45181250

VITA VM 13 ENAMEL, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
ENL	B4519112	B4519150	B45191250
END	B4519212	B4519250	B45192250

VITA VM 13 NEUTRAL, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
NT	B4520112	B4520150	B45201250

VITA VM 13 CORRECTIVE, 12 g

Description	Prod. no. 12 g
COR1	B4522112
COR2	B4522212
COR3	B4522312

VITA VM 13 GINGIVA, 12 g

Description	Prod. no. 12 g
G1	B4523112
G2	B4523212
G3	B4523312
G4	B4523412
G5	B4523512

VITA VM 13 CHROMA PLUS, 12 g

Description	Prod. no. 12 g
CP1	B4524112
CP2	B4524212
CP3	B4524312
CP4	B4524412
CP5	B4524512

VITA VM 13 WASH OPAQUE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
W0	B4525912	B4525950

VITA VM 13 MARGIN, 12 g

Description	Prod. no. 12 g
M1	B4526112
M4	B4526412
M5	B4526512
M7	B4526712
M8	B4526812
M9	B4526912

VITA VM 13 GINGIVA OPAQUE, 12 g

Description	Prod. no. 12 g
GOL	B4527112
GOD	B4527312

VITA VM 13 COLOR OPAQUE, 12 g

Description	Prod. no. 12 g
C01	B4528112
C02	B4528212
C03	B4528312

VITA VM 13 SUN OPAQUE, 12 g

Description	Prod. no. 12 g
S01	B4529112
S02	B4529212
S03	B4529312

VITA VM 13 WASH OPAQUE PASTE, 7 g

Description	Prod. no. 7 g
W0	B453597

VITA VM 13 SUN DENTINE, 12 g

Description	Prod. no. 12 g
SD1	B4536112
SD2	B4536212
SD3	B4536312

VITA VM 13 GINGIVA OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
G0L	B453715
G0D	B453735

VITA VM 13 COLOR OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
C01	B453815
C02	B453825
C03	B453835

VITA VM 13 SUN OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
S01	B453915
S02	B453925
S03	B453935

VITA VMK Master®

For full veneering of metal frameworks in the conventional CTE range



Product description

VITA VMK Master is a natural feldspar veneering ceramic for conventional bonding alloys (approx. 13.8 – 15.2 $\mu\text{m}/\text{mK}$).

Regardless of whether the conventional alloys are cast, milled or sintered.

Precise and reliable processing.

A master in terms of reliability, suitability for everyday use and esthetics – that's VITA VMK Master.

Features/Benefits

- A material that makes no compromises when it comes to reliability and esthetics
- Precise and reliable processing, allowing it to be learned quickly
- Excellent processing characteristics, such as improved modeling properties, simple handling and minimum shrinkage behavior
- Conventional, established dentin-incisal layering system
- Comprehensive range of additional materials for individualization
- Highly esthetic restorations thanks to brilliant colors
- Available in VITA SYSTEM 3D-Master and VITA classical A1 – D4 shades

VITA VMK Master® kits

VITA VMK Master® VITA classical A1-D4®



VITA VMK Master® ONE COLOR KIT A2

One-shade kit for initial tests

BVMKOCSA2V2

Quantity	Content	Material
1	12 g	VITA VMK Master OPAQUE A2
1	12 g	VITA VMK Master OPAQUE DENTINE A2
1	12 g	VITA VMK Master DENTINE A2
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1	20 ml	VITA AKZENT Plus POWDER FLUID

VITA VMK Master® STANDARD SET POWDER VITA classical A1-D4®

Basic set for the reproduction of the VITA classical A1-D4 shades

BVMKSETCV1

Quantity	Content	Material
1	12 g	VITA VMK Master WASH OPAQUE WO
16	12 g	VITA VMK Master OPAQUE A1 - D4
16	12 g	VITA VMK Master OPAQUE DENTINE A1 - D4
16	12 g	VITA VMK Master DENTINE A1 - D4
2	12 g	VITA VMK Master ENAMEL EN1 - EN2
1	12 g	VITA VMK Master TRANSLUCENT T4
1	12 g	VITA VMK Master WINDOW WIN
3	12 g	VITA VMK Master CERVICAL CE1 - CE3
5	12 g	VITA VMK Master NECK N1 - N5
2	12 g	VITA VMK Master CORRECTIVE COR1 - COR2
1	50 ml	VITA MODELLING FLUID
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VMK MASTER classical
1		VITA classical A1-D4 shade guide
1		Shade guide for VITA VMK MASTER STANDARD

VITA VMK Master® STANDARD SET PASTE VITA classical A1-D4®
Basic set for the reproduction of the VITA classical A1-D4 shades
BVMKSETPCV1

Quantity	Content	Material
1	7 g	VITA VMK Master WASH OPAQUE PASTE W0
16	5 g	VITA VMK Master OPAQUE PASTE A1 - D4
16	12 g	VITA VMK Master OPAQUE DENTINE A1 - D4
16	12 g	VITA VMK Master DENTINE A1 - D4
2	12 g	VITA VMK Master ENAMEL EN1 - EN2
1	12 g	VITA VMK Master TRANSLUCENT T4
1	12 g	VITA VMK Master WINDOW WIN
3	12 g	VITA VMK Master CERVICAL CE1 - CE3
5	12 g	VITA VMK Master NECK N1 - N5
2	12 g	VITA VMK Master CORRECTIVE COR1 - COR2
1	50 ml	VITA MODELLING FLUID
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA PASTE OPAQUE LIQUID
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VMK MASTER classical
1		VITA classical A1-D4 shade guide
1		Shade guide for VITA VMK MASTER STANDARD



VITA VMK Master® STARTER SET classical
Two-color set
BVMKSSC

Quantity	Content	Material
2	12 g	VITA VMK Master OPAQUE A2 - A3
2	12 g	VITA VMK Master OPAQUE DENTINE A2 - A3
2	12 g	VITA VMK Master DENTINE A2 - A3
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1		VITA shade indicator stick A2
1		VITA shade indicator stick A3

VITA VMK Master® STARTER SET PASTE VITA classical A1-D4®

Two-color set

BVMKSSCP

Quantity	Content	Material
2	5 g	VITA VMK Master OPAQUE PASTE A2 - A3
2	12 g	VITA VMK Master OPAQUE DENTINE A2 - A3
2	12 g	VITA VMK Master DENTINE A2 - A3
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA PASTE OPAQUE LIQUID
1		VITA shade indicator stick A2
1		VITA shade indicator stick A3
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		Flat brush, for PASTE OPAQUE

VITA VMK Master® VITA SYSTEM 3D-MASTER®



VITA VMK Master® 10-COLOR-SET POWDER VITA SYSTEM 3D-MASTER® Reduced basic set with the ten most common 3D-MASTER shades BVMK10S3D

Quantity	Content	Material
1	12 g	VITA VMK Master WASH OPAQUE WO
4	12 g	VITA VMK Master OPAQUE OP1 - OP4
10	12 g	VITA VMK Master OPAQUE DENTINE 1M2, 2L1.5, 2M1 - 2M2, 2R2.5, 3L1.5, 3M1 - 3M3, 4M2
10	12 g	VITA VMK Master DENTINE 1M2, 2L1.5, 2M1 - 2M2, 2R2.5, 3L1.5, 3M1 - 3M3, 4M2
2	12 g	VITA VMK Master ENAMEL EN1, EN3
1	12 g	VITA VMK Master TRANSLUCENT T4
1	12 g	VITA VMK Master WINDOW WIN
3	12 g	VITA VMK Master CERVICAL CE1 - CE3
5	12 g	VITA VMK Master NECK N1 - N5
2	12 g	VITA VMK Master CORRECTIVE COR1 - COR2
1	50 ml	VITA MODELLING FLUID
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VMK MASTER 3D-MASTER
1		VITA Linearguide 3D-MASTER
1		Shade guide for VITA VMK MASTER STANDARD

VITA VMK Master® 10-COLOR-SET PASTE VITA SYSTEM 3D-MASTER®
Reduced basic set with the ten most common 3D-MASTER shades
BVMK10SP3D

Quantity	Content	Material
1	7 g	VITA VMK Master WASH OPAQUE PASTE W0
4	5 g	VITA VMK Master OPAQUE PASTE OP1 - OP4
10	12 g	VITA VMK Master OPAQUE DENTINE 1M2, 2L1.5, 2M1 - 2M2, 2R2.5, 3L1.5, 3M1 - 3M3, 4M2
10	12 g	VITA VMK Master DENTINE 1M2, 2L1.5, 2M1 - 2M2, 2R2.5, 3L1.5, 3M1 - 3M3, 4M2
2	12 g	VITA VMK Master ENAMEL EN1, EN3
1	12 g	VITA VMK Master TRANSLUCENT T4
1	12 g	VITA VMK Master WINDOW WIN
3	12 g	VITA VMK Master CERVICAL CE1 - CE3
5	12 g	VITA VMK Master NECK N1 - N5
2	12 g	VITA VMK Master CORRECTIVE COR1 - COR2
1	50 ml	VITA MODELLING FLUID
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA PASTE OPAQUE LIQUID
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VMK MASTER 3D-MASTER
1		VITA Linearguide 3D-MASTER
1		Shade guide for VITA VMK MASTER STANDARD



VITA VMK Master® BLEACHED COLOR SET POWDER VITA SYSTEM
3D-MASTER®
Very light shades to reproduce bleached teeth
BVMKBCS

Quantity	Content	Material
1	12 g	VITA VMK Master OPAQUE OP0
3	12 g	VITA VMK Master DENTINE OM1 - OM3
3	12 g	VITA VMK Master OPAQUE DENTINE OM1 - OM3
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1		VITA BLEACHED SHADE GUIDE, Shade Group 0

VITA VMK Master® BLEACHED COLOR SET PASTE VITA SYSTEM 3D-MASTER®

Very light shades to reproduce bleached teeth

BVMKBCSP

Quantity	Content	Material
1	5 g	VITA VMK Master OPAQUE PASTE OP0
3	12 g	VITA VMK Master DENTINE 0M1 - 0M3
3	12 g	VITA VMK Master OPAQUE DENTINE 0M1 - 0M3
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA PASTE OPAQUE LIQUID
1		VITA BLEACHED SHADE GUIDE, Shade Group 0
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		Flat brush, for PASTE OPAQUE



VITA VMK Master® ONE COLOR KIT 2M2

One-shade kit for initial tests

BVMKOC2M2V2

Quantity	Content	Material
1	12 g	VITA VMK Master OPAQUE OP2
1	12 g	VITA VMK Master OPAQUE DENTINE 2M2
1	12 g	VITA VMK Master DENTINE 2M2
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1	20 ml	VITA AKZENT Plus POWDER FLUID

VITA VMK Master® STANDARD SET POWDER VITA SYSTEM

3D-MASTER®

Basic set for the reproduction of the 3D-MASTER shades

BVMKSET3DV1

Quantity	Content	Material
1	12 g	VITA VMK Master WASH OPAQUE WO
5	12 g	VITA VMK Master OPAQUE OP1 - OP5
26	12 g	VITA VMK Master DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
26	12 g	VITA VMK Master OPAQUE DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
2	12 g	VITA VMK Master ENAMEL EN1, EN3
1	12 g	VITA VMK Master TRANSLUCENT T4
1	12 g	VITA VMK Master WINDOW WIN
3	12 g	VITA VMK Master CERVICAL CE1 - CE3
5	12 g	VITA VMK Master NECK N1 - N5
2	12 g	VITA VMK Master CORRECTIVE COR1 - COR2
1	50 ml	VITA MODELLING FLUID
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VMK MASTER 3D-MASTER
1		VITA Linearguide 3D-MASTER
1		Shade guide for VITA VMK MASTER STANDARD

VITA VMK Master® STANDARD SET PASTE VITA SYSTEM 3D-MASTER®
Basic set for the reproduction of the 3D-MASTER shades
BVMKSETP3DV1

Quantity	Content	Material
1	7 g	VITA VMK Master WASH OPAQUE PASTE WO
5	5 g	VITA VMK Master OPAQUE PASTE OP1 - OP5
26	12 g	VITA VMK Master DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
26	12 g	VITA VMK Master OPAQUE DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
2	12 g	VITA VMK Master ENAMEL EN1, EN3
1	12 g	VITA VMK Master TRANSLUCENT T4
1	12 g	VITA VMK Master WINDOW WIN
3	12 g	VITA VMK Master CERVICAL CE1 - CE3
5	12 g	VITA VMK Master NECK N1 - N5
2	12 g	VITA VMK Master CORRECTIVE COR1 - COR2
1	50 ml	VITA MODELLING FLUID
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA PASTE OPAQUE LIQUID
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		VITA Modisol®
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VMK MASTER 3D-MASTER
1		VITA Linearguide 3D-MASTER
1		Shade guide for VITA VMK MASTER STANDARD



VITA VMK Master® STARTER SET POWDER VITA SYSTEM 3D-MASTER®
Two-color set
BVMKSS3D

Quantity	Content	Material
2	12 g	VITA VMK Master OPAQUE OP2 - OP3
2	12 g	VITA VMK Master DENTINE 2M2, 3M2
2	12 g	VITA VMK Master OPAQUE DENTINE 2M2, 3M2
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	50 ml	VITA MODELLING FLUID RS
1	50 ml	VITA OPAQUE FLUID
1		VITA shade indicator stick 2M2
1		VITA shade indicator stick 3M2

VITA VMK Master® STARTER SET PASTE VITA SYSTEM 3D-MASTER®
Two-color set
BVMKSP3D

Quantity	Content	Material
2	5 g	VITA VMK Master OPAQUE PASTE OP2 - OP3
2	12 g	VITA VMK Master DENTINE 2M2, 3M2
2	12 g	VITA VMK Master OPAQUE DENTINE 2M2, 3M2
1	12 g	VITA VMK Master ENAMEL EN1
1	12 g	VITA VMK Master TRANSLUCENT T4
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA PASTE OPAQUE LIQUID
1		VITA shade indicator stick 2M2
1		VITA shade indicator stick 3M2
1		Flat brush, for PASTE OPAQUE
1		Mixing spatula, Acrylic for PASTE OPAQUE

VITA VMK Master® Zusatzmassen
compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades



VITA VMK Master® ADDITIONAL SET
Additional materials for achieving special shade effects and individualizing the restoration
BVMKAS

Quantity	Content	Material
8	12 g	VITA VMK Master TRANSLUCENT T1 - T8
8	12 g	VITA VMK Master DENTINE MODIFIER DM1 - DM8
6	12 g	VITA VMK Master LUMINARY LM1 - LM6
6	12 g	VITA VMK Master MARGIN MN, M1 - M5
1	12 g	VITA VMK Master PEARL TRANSLUCENT PLT1
1	12 g	VITA VMK Master OPAL TRANSLUCENT OT1
1	50 ml	VITA MODELLING FLUID
1		Shade guide for VITA VMK MASTER TRANSLUCENT
1		Shade guide for VITA VMK MASTER MARGIN / GINGIVA
1		Shade guide for VITA VMK MASTER DENTINE MODIFIER / LUMINARY



VITA VMK Master® GINGIVA SET POWDER
Natural-looking gingiva pastes
BVMKGS

Quantity	Content	Material
5	12 g	VITA VMK Master GINGIVA G1 - G5
2	12 g	VITA VMK Master GINGIVA OPAQUE GOL, GOD
1		Shade guide for VITA VMK MASTER MARGIN / GINGIVA

VITA VMK Master® GINGIVA SET PASTE
Natural-looking gingiva pastes
BVMKGSP

Quantity	Content	Material
5	12 g	VITA VMK Master GINGIVA G1 - G5
2	5 g	VITA VMK Master GINGIVA OPAQUE PASTE GOL, GOD
1		Shade guide for VITA VMK MASTER MARGIN / GINGIVA

VITA VMK Master® products

VITA VMK Master® VITA classical A1-D4®

VITA VMK Master OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
A1	B483015
A2	B483025
A3	B483035
A3.5	B483045
A4	B483055
B1	B483065
B2	B483075
B3	B483085
B4	B483095
C1	B483105
C2	B483115
C3	B483125
C4	B483135
D2	B483145
D3	B483155
D4	B483165

VITA VMK Master OPAQUE, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
A1	B4840112	B4840150	B48401250
A2	B4840212	B4840250	B48402250
A3	B4840312	B4840350	B48403250
A3.5	B4840412	B4840450	-
A4	B4840512	B4840550	-
B1	B4840612	B4840650	-
B2	B4840712	B4840750	-
B3	B4840812	B4840850	-
B4	B4840912	B4840950	-
C1	B4841012	B4841050	-
C2	B4841112	B4841150	-
C3	B4841212	B4841250	-
C4	B4841312	B4841350	-
D2	B4841412	B4841450	-
D3	B4841512	B4841550	-
D4	B4841612	B4841650	-

VITA VMK Master OPAQUE DENTINE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
A1	B4843112	B4843150
A2	B4843212	B4843250
A3	B4843312	B4843350
A3.5	B4843412	B4843450
A4	B4843512	B4843550
B1	B4843612	B4843650
B2	B4843712	B4843750
B3	B4843812	B4843850
B4	B4843912	B4843950
C1	B4844012	B4844050
C2	B4844112	B4844150
C3	B4844212	B4844250
C4	B4844312	B4844350
D2	B4844412	B4844450
D3	B4844512	B4844550
D4	B4844612	B4844650

VITA VMK Master DENTINE, 12 g / 50 g / 250 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
A1	EN1	B4845112	B4845150	B48451250
A2	EN1	B4845212	B4845250	B48452250
A3	EN1	B4845312	B4845350	B48453250
A3.5	EN2	B4845412	B4845450	B48454250
A4	EN2	B4845512	B4845550	-
B1	EN2	B4845612	B4845650	B48456250
B2	EN2	B4845712	B4845750	B48457250
B3	EN2	B4845812	B4845850	-
B4	EN2	B4845912	B4845950	-
C1	EN2	B4846012	B4846050	B48460250
C2	EN2	B4846112	B4846150	B48461250
C3	EN1	B4846212	B4846250	B48462250
C4	EN1	B4846312	B4846350	-
D2	EN2	B4846412	B4846450	B48464250
D3	EN2	B4846512	B4846550	-
D4	EN2	B4846612	B4846650	-

VITA VMK Master® VITA SYSTEM 3D-MASTER®

VITA VMK Master DENTINE, 12 g / 50 g / 250 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
0M1	EN1	B4803112	B4803150	-
0M2	EN1	B4803212	B4803250	-
0M3	EN1	B4803312	B4803350	-
1M1	EN1	B4803412	B4803450	B48034250
1M2	EN1	B4803512	B4803550	B48035250
2L1.5	EN1	B4803612	B4803650	B48036250
2L2.5	EN1	B4803712	B4803750	-
2M1	EN1	B4803812	B4803850	B48038250
2M2	EN1	B4803912	B4803950	B48039250
2M3	EN1	B4804012	B4804050	-
2R1.5	EN1	B4804112	B4804150	-
2R2.5	EN1	B4804212	B4804250	B48042250
3L1.5	EN1	B4804312	B4804350	-
3L2.5	EN1	B4804412	B4804450	-
3M1	EN1	B4804512	B4804550	-
3M2	EN1	B4804612	B4804650	B48046250
3M3	EN1	B4804712	B4804750	-

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
3R1.5	EN1	B4804812	B4804850	-
3R2.5	EN3	B4804912	B4804950	-
4L1.5	EN1	B4805012	B4805050	-
4L2.5	EN1	B4805112	B4805150	-
4M1	EN1	B4805212	B4805250	-
4M2	EN3	B4805312	B4805350	-
4M3	EN3	B4805412	B4805450	-
4R1.5	EN1	B4805512	B4805550	-
4R2.5	EN3	B4805612	B4805650	-
5M1	EN1	B4805712	B4805750	-
5M2	EN3	B4805812	B4805850	-
5M3	EN3	B4805912	B4805950	-

VITA VMK Master OPAQUE DENTINE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
0M1	B4806112	B4806150
0M2	B4806212	B4806250
0M3	B4806312	B4806350
1M1	B4806412	B4806450
1M2	B4806512	B4806550
2L1.5	B4806612	B4806650
2L2.5	B4806712	B4806750
2M1	B4806812	B4806850
2M2	B4806912	B4806950
2M3	B4807012	B4807050
2R1.5	B4807112	B4807150
2R2.5	B4807212	B4807250
3L1.5	B4807312	B4807350
3L2.5	B4807412	B4807450
3M1	B4807512	B4807550
3M2	B4807612	B4807650
3M3	B4807712	B4807750
3R1.5	B4807812	B4807850
3R2.5	B4807912	B4807950
4L1.5	B4808012	B4808050
4L2.5	B4808112	B4808150
4M1	B4808212	B4808250
4M2	B4808312	B4808350
4M3	B4808412	B4808450

Description	Prod. no. 12 g	Prod. no. 50 g
4R1.5	B4808512	B4808550
4R2.5	B4808612	B4808650
5M1	B4808712	B4808750
5M2	B4808812	B4808850
5M3	B4808912	B4808950

VITA VMK Master OPAQUE, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
OP0	B4825012	B4825050	-
OP1	B4825112	B4825150	B48251250
OP2	B4825212	B4825250	B48252250
OP3	B4825312	B4825350	B48253250
OP4	B4825412	B4825450	-
OP5	B4825512	B4825550	-

VITA VMK Master OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
OP0	B483505
OP1	B483515
OP2	B483525
OP3	B483535
OP4	B483545
OP5	B483555

VITA VMK Master® Zusatzmassen
compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades
VITA VMK Master TRANSLUCENT, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
T1	B4809112	-	-
T2	B4809212	-	-
T3	B4809312	-	-
T4	B4809412	B4809450	B48094250
T5	B4809512	-	-
T6	B4809612	-	-
T7	B4809712	-	-
T8	B4809812	-	-

VITA VMK Master DENTINE MODIFIER, 12 g

Description	Prod. no. 12 g
DM1	B4812112
DM2	B4812212
DM3	B4812312
DM4	B4812412
DM5	B4812512
DM6	B4812612
DM7	B4812712
DM8	B4812812

VITA VMK Master LUMINARY, 12 g

Description	Prod. no. 12 g
LM1	B4814112
LM2	B4814212
LM3	B4814312
LM4	B4814412
LM5	B4814512
LM6	B4814612

VITA VMK Master PEARL TRANSLUCENT, 12 g

Description	Prod. no. 12 g
PLT1	B4815112

VITA VMK Master OPAL TRANSLUCENT, 12 g

Description	Prod. no. 12 g
OT1	B4817112

VITA VMK Master WINDOW, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
WIN	B4818112	B4818150	B48181250

VITA VMK Master ENAMEL, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
EN1	B4819112	B4819150	B48191250
EN2	B4819212	B4819250	B48192250
EN3	B4819312	B4819350	B48193250

VITA VMK Master CERVICAL, 12 g

Description	Prod. no. 12 g
CE1	B4820112
CE2	B4820212
CE3	B4820312

VITA VMK Master CORRECTIVE, 12 g

Description	Prod. no. 12 g
COR1	B4822112
COR2	B4822212

VITA VMK Master GINGIVA, 12 g

Description	Prod. no. 12 g
G1	B4823112
G2	B4823212
G3	B4823312
G4	B4823412
G5	B4823512

VITA VMK Master NECK, 12 g

Description	Prod. no. 12 g
N1	B4824112
N2	B4824212
N3	B4824312
N4	B4824412
N5	B4824512

VITA VMK Master WASH OPAQUE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
W0	B4825912	B4825950

VITA VMK Master MARGIN, 12 g

Description	Prod. no. 12 g
MN	B4826012
M1	B4826112
M2	B4826212
M3	B4826312
M4	B4826412
M5	B4826512

VITA VMK Master GINGIVA OPAQUE, 12 g

Description	Prod. no. 12 g
GOL	B4827112
GOD	B4827312

VITA VMK Master WASH OPAQUE PASTE, 7 g

Description	Prod. no. 7 g
WO	B483597

VITA VMK Master GINGIVA OPAQUE PASTE, 5 g

Description	Prod. no. 5 g
GOL	B483715
GOD	B483735

VITA NP BOND

Effective CTE buffer between the ceramic and the metal alloy



Product description

Bonder for greater reliability and improved esthetics when veneering non-precious metal dental alloys. Through reliable bonding between the metal framework and the veneering ceramic (such as VITA VMK MASTER and VITA VM 13), VITA NP BOND helps prevent errors occurring as a result of CTE tension between the alloy and the ceramic.

VITA NP BOND (NP: non-precious) improves bonding between metal frameworks and veneering ceramics, particularly in the case of non-precious metal dental alloys, thus providing dental technicians with greater reliability in fabricating veneers that offer superior long-term stability and esthetics.

Features/Benefits

- Effective CTE buffer between the ceramic and the metal alloy, particularly in the case of non-precious metal alloys, provides for a stable bond between the framework and the ceramic
- Improves shade reproduction
- Saves time as firing with WASH OPAQUE is no longer necessary
- Thanks to the improved surface structure following BONDER firing, the opaquer can be applied more easily, effectively and evenly

VITA NP BOND products

VITA NP BOND

Description	Prod. no.
VITA NP BOND PASTE, 5 g	BNPB5

VITA LUMEX® AC

Ideal shade fidelity. Excellent light dynamics. Precise processing.



Product description

- VITA LUMEX AC is a leucite-reinforced, glass-ceramic veneering system with an innovative material formula for vivid and reliable results with shade accuracy
- VITA LUMEX AC includes:
- GINGIVA, OPAQUE, OPAQUE DENTINE, DENTINE and ENAMEL materials, as well as a multitude of effect materials (e.g., OPAL TRANSLUCENT, FLUO INTENSE and much more)
- ideally matched system components such as brushes, firing trays, glazing materials, etc.

Features/Benefits

- Reliably reproduces tooth shades, as well as a natural play of shade and light, thanks to ceramic materials with very good shade fidelity and excellent light dynamics
- Process veneering ceramic efficiently and precisely, thanks to the ideal modeling properties, outstanding shape stability and universal applicability for all substructure ceramics
- Produce durable restorations, thanks to the high fracture resistance of the leucite-enhanced glass ceramic matrix and a reliable bond to the substructure

Not available in all markets. Official approval pending in some markets.

Indications

For the veneering of all common full ceramic substructure materials (zirconia, LiSi-glass and feldspar ceramics) and for the production of reconstructions without a substructure (e.g., veneers)

VITA LUMEX® AC kits

VITA LUMEX® AC VITA classical A1-D4®

VITA LUMEX® AC STANDARD KIT VITA classical A1-D4®

For standard veneering

BLSKCV1

Quantity	Content	Material
16	12 g	VITA LUMEX AC OPAQUE DENTINE A1 - D4
16	12 g	VITA LUMEX AC DENTINE A1 - D4
5	12 g	VITA LUMEX AC ENAMEL light, medium, intense, clear, fog
2	12 g	VITA LUMEX AC CORRECTIVE neutral, desert
1	50 ml	VITA LUMEX AC MODELLING LIQUID
1	12 g	VITA LUMEX AC MARGIN, straw-yellow
1	12 g	VITA LUMEX AC MARGIN, corn-yellow
1	12 g	VITA LUMEX AC TRANSLUCENT, light-blonde
1	5 g	VITA AKZENT Plus GLAZE LT POWDER
1	20 ml	VITA AKZENT Plus POWDER FLUID
1		Paint brush no. 3/0 synthetic
1		Coating brush no.8 synthetic
1		Firing tray G (pack of 2 trays and 20 pins)

VITA LUMEX® AC STARTER KIT VITA classical A1-D4®

For the first introduction

BLSTKCV1

Quantity	Content	Material
4	12 g	VITA LUMEX AC OPAQUE DENTINE A1 - A3.5
4	12 g	VITA LUMEX AC DENTINE A1 - A3.5
4	12 g	VITA LUMEX AC ENAMEL light, medium, clear, fog
2	12 g	VITA LUMEX AC CORRECTIVE neutral, desert
1	12 g	VITA LUMEX AC PEARL shell
1	12 g	VITA LUMEX AC TRANSLUCENT, smoky-white
1	12 g	VITA LUMEX AC TRANSLUCENT, light-blond
1	12 g	VITA LUMEX AC TRANSLUCENT waterdrop
1	12 g	VITA LUMEX AC TRANSLUCENT, foggy-grey
1	12 g	VITA LUMEX AC FLUO INTENSE, arctic-white
2	12 g	VITA LUMEX AC FLUO INTENSE cream, cappuccino
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-neutral
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-sky
1	50 ml	VITA LUMEX AC MODELLING LIQUID
1	5 g	VITA AKZENT Plus GLAZE LT POWDER
1	20 ml	VITA AKZENT Plus POWDER FLUID
1		Paint brush no. 3/0 synthetic
1		Coating brush no.8 synthetic
1		Firing tray G (pack of 2 trays and 20 pins)

VITA LUMEX® AC additional powder



VITA LUMEX® AC CUT-BACK KIT For cut-back and micro-veneering BLCBK

Quantity	Content	Material
5	12 g	VITA LUMEX AC ENAMEL light, medium, intense, clear, fog
1	12 g	VITA LUMEX AC TRANSLUCENT, smoky-white
1	12 g	VITA LUMEX AC TRANSLUCENT, light-blonde
1	12 g	VITA LUMEX AC TRANSLUCENT, misty-rose
2	12 g	VITA LUMEX AC TRANSLUCENT sunlight, waterdrop
1	12 g	VITA LUMEX AC TRANSLUCENT, sun-intense
1	12 g	VITA LUMEX AC TRANSLUCENT, deep-blue
1	12 g	VITA LUMEX AC TRANSLUCENT, foggy-grey
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-neutral
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-sky
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-azure
2	12 g	VITA LUMEX AC CORRECTIVE neutral, desert
1	50 ml	VITA LUMEX AC MODELLING LIQUID
1	5 g	VITA AKZENT Plus GLAZE LT POWDER
1	20 ml	VITA AKZENT Plus POWDER FLUID



VITA LUMEX® AC EXPERT KIT

**For individual veneering
BLEK**

Quantity	Content	Material
1	12 g	VITA LUMEX AC TRANSLUCENT, smoky-white
1	12 g	VITA LUMEX AC TRANSLUCENT, light-blonde
1	12 g	VITA LUMEX AC TRANSLUCENT, misty-rose
2	12 g	VITA LUMEX AC TRANSLUCENT sunlight, waterdrop
1	12 g	VITA LUMEX AC TRANSLUCENT, sun-intense
1	12 g	VITA LUMEX AC TRANSLUCENT, deep-blue
1	12 g	VITA LUMEX AC TRANSLUCENT, foggy-grey
3	12 g	VITA LUMEX AC CHROMA INTENSE ivory, almond, hazelnut
1	12 g	VITA LUMEX AC FLUO INTENSE, arctic-white
4	12 g	VITA LUMEX AC FLUO INTENSE cream, cappuccino, sand, sesame
1	12 g	VITA LUMEX AC DENTINE MODIFIER, cloudy-white
4	12 g	VITA LUMEX AC DENTINE MODIFIER caramel, honey, copper, brown
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-neutral
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-sky
1	12 g	VITA LUMEX AC OPAL TRANSLUCENT, opal-azure
1	12 g	VITA LUMEX AC MAMELON saffron
1	12 g	VITA LUMEX AC MAMELON, honey-melon
1	12 g	VITA LUMEX AC PEARL shell
1	50 ml	VITA LUMEX AC MODELLING LIQUID
1	5 g	VITA AKZENT Plus GLAZE LT POWDER
1	20 ml	VITA AKZENT Plus POWDER FLUID



VITA LUMEX® AC GINGIVA KIT

**For the reconstruction of gingival areas
BLGK**

Quantity	Content	Material
1	12 g	VITA LUMEX AC GINGIVA, pale-papilla
1	12 g	VITA LUMEX AC GINGIVA, light-rose
4	12 g	VITA LUMEX AC GINGIVA nectarine, grapefruit, rosewood, purple
1	12 g	VITA LUMEX AC GINGIVA, deep-red
1	12 g	VITA LUMEX AC GINGIVA, dark-red
1	5 g	VITA AKZENT Plus GLAZE LT POWDER
1	20 ml	VITA AKZENT Plus POWDER FLUID

VITA LUMEX® AC VITA SYSTEM 3D-MASTER®

VITA LUMEX® AC STANDARD KIT VITA SYSTEM 3D-MASTER®

For standard veneering

BLSKV1

Quantity	Content	Material
29	12 g	VITA LUMEX AC OPAQUE DENTINE 0M1 - 0M3, 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
29	12 g	VITA LUMEX AC DENTINE 0M1 - 0M3, 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
5	12 g	VITA LUMEX AC ENAMEL light, medium, intense, clear, fog
2	12 g	VITA LUMEX AC CORRECTIVE neutral, desert
1	12 g	VITA LUMEX AC MARGIN, straw-yellow
1	12 g	VITA LUMEX AC MARGIN, corn-yellow
1	12 g	VITA LUMEX AC TRANSLUCENT, light-blond
1	50 ml	VITA LUMEX AC MODELLING LIQUID
1	5 g	VITA AKZENT Plus GLAZE LT POWDER
1	20 ml	VITA AKZENT Plus POWDER FLUID
1		Paint brush no. 3/0 synthetic
1		Coating brush no.8 synthetic
1		Firing tray G (pack of 2 trays and 20 pins)

VITA LUMEX® AC products

VITA LUMEX® AC VITA classical A1-D4®

VITA LUMEX AC OPAQUE DENTINE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
A1	B5343112	B5343150
A2	B5343212	B5343250
A3	B5343312	B5343350
A3.5	B5343412	B5343450
A4	B5343512	B5343550
B1	B5343612	B5343650
B2	B5343712	B5343750
B3	B5343812	B5343850
B4	B5343912	B5343950
C1	B5344012	B5344050
C2	B5344112	B5344150
C3	B5344212	B5344250
C4	B5344312	B5344350
D2	B5344412	B5344450
D3	B5344512	B5344550
D4	B5344612	B5344650

VITA LUMEX AC DENTINE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
A1	B5345112	B5345150
A2	B5345212	B5345250
A3	B5345312	B5345350
A3.5	B5345412	B5345450
A4	B5345512	B5345550
B1	B5345612	B5345650
B2	B5345712	B5345750
B3	B5345812	B5345850
B4	B5345912	B5345950
C1	B5346012	B5346050
C2	B5346112	B5346150
C3	B5346212	B5346250
C4	B5346312	B5346350
D2	B5346412	B5346450
D3	B5346512	B5346550
D4	B5346612	B5346650

VITA LUMEX® AC additional powder

VITA LUMEX AC TRANSLUCENT, 12 g

Description	Prod. no. 12 g
smoky-white	B5309112
light-blonde	B5309212
misty-rose	B5309312
sunlight	B5309412
sun-intense	B5309512
waterdrop	B5309612
deep-blue	B5309712
foggy-grey	B5309812

VITA LUMEX AC DENTINE MODIFIER, 12 g

Description	Prod. no. 12 g
cloudy-white	B5312112
caramel	B5312212
honey	B5312312
copper	B5312412
brown	B5312512

VITA LUMEX AC FLUO INTENSE, 12 g

Description	Prod. no. 12 g
arctic-white	B5314112
cream	B5314212
cappuccino	B5314312
sand	B5314412
sesame	B5314512

VITA LUMEX AC PEARL, 12 g

Description	Prod. no. 12 g
shell	B5315112

VITA LUMEX AC MAMELON, 12 g

Description	Prod. no. 12 g
saffron	B5316112
honey-melon	B5316212

VITA LUMEX AC OPAL TRANSLUCENT, 12 g

Description	Prod. no. 12 g
opal-neutral	B5317112
opal-sky	B5317212
opal-azure	B5317312

VITA LUMEX AC ENAMEL, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
clear	B5318112	B5318150
light	B5319112	B5319150
medium	B5319212	B5319250
intense	B5319312	B5319350
fog	B5320112	B5320150

VITA LUMEX AC CORRECTIVE, 12 g

Description	Prod. no. 12 g
neutral	B5322112
desert	B5322212

VITA LUMEX AC GINGIVA, 12 g

Description	Prod. no. 12 g
pale-papilla	B5323112
light-rose	B5323212
nectarine	B5323312
grapefruit	B5323412
rosewood	B5323512
purple	B5323612
deep-red	B5323712
dark-red	B5323812

VITA LUMEX AC MARGIN, 12 g

Description	Prod. no. 12 g
straw-yellow	B5326112
corn-yellow	B5326212

VITA LUMEX AC CHROMA INTENSE, 12 g

Description	Prod. no. 12 g
ivory	B5351112
almond	B5351212
hazelnut	B5351312

VITA LUMEX® AC VITA SYSTEM 3D-MASTER®

VITA LUMEX AC DENTINE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
0M1	B5303112	B5303150
0M2	B5303212	B5303250
0M3	B5303312	B5303350
1M1	B5303412	B5303450
1M2	B5303512	B5303550
2L1.5	B5303612	B5303650
2L2.5	B5303712	B5303750
2M1	B5303812	B5303850
2M2	B5303912	B5303950
2M3	B5304012	B5304050
2R1.5	B5304112	B5304150
2R2.5	B5304212	B5304250
3L1.5	B5304312	B5304350
3L2.5	B5304412	B5304450

Description	Prod. no. 12 g	Prod. no. 50 g
3M1	B5304512	B5304550
3M2	B5304612	B5304650
3M3	B5304712	B5304750
3R1.5	B5304812	B5304850
3R2.5	B5304912	B5304950
4L1.5	B5305012	B5305050
4L2.5	B5305112	B5305150
4M1	B5305212	B5305250
4M2	B5305312	B5305350
4M3	B5305412	B5305450
4R1.5	B5305512	B5305550
4R2.5	B5305612	B5305650
5M1	B5305712	B5305750
5M2	B5305812	B5305850
5M3	B5305912	B5305950

VITA LUMEX AC OPAQUE DENTINE, 12 g / 50 g

Description	Prod. no. 12 g	Prod. no. 50 g
0M1	B5306112	B5306150
0M2	B5306212	B5306250
0M3	B5306312	B5306350
1M1	B5306412	B5306450
1M2	B5306512	B5306550
2L1.5	B5306612	B5306650
2L2.5	B5306712	B5306750
2M1	B5306812	B5306850
2M2	B5306912	B5306950
2M3	B5307012	B5307050
2R1.5	B5307112	B5307150
2R2.5	B5307212	B5307250
3L1.5	B5307312	B5307350
3L2.5	B5307412	B5307450
3M1	B5307512	B5307550
3M2	B5307612	B5307650
3M3	B5307712	B5307750
3R1.5	B5307812	B5307850
3R2.5	B5307912	B5307950
4L1.5	B5308012	B5308050
4L2.5	B5308112	B5308150

Description	Prod. no. 12 g	Prod. no. 50 g
4M1	B5308212	B5308250
4M2	B5308312	B5308350
4M3	B5308412	B5308450
4R1.5	B5308512	B5308550
4R2.5	B5308612	B5308650
5M1	B5308712	B5308750
5M2	B5308812	B5308850
5M3	B5308912	B5308950

VITA LUMEX AC OPAQUE, 12 g

Description	Prod. no. 12 g
opaque-0	B5325012
opaque-1	B5325112
opaque-2	B5325212
opaque-3	B5325312
opaque-4	B5325412
opaque-5	B5325512

VITAVM®11

Fine-structure ceramic for individualization of VITA SUPRINITY PC*



Product description

VITA VM 11 is a low-melting, fine-structure feldspar ceramic that is adapted specifically for the individualization of restorations made of zirconia-reinforced lithium silicate glass ceramic. The optimum correspondence of both CTE ranges guarantees a bond that is tension-free, and in turn also reliable.

In combination with lithium silicate-type material, you can begin layering directly without liner firing or washbake, thanks to the excellent surface wettability of the new glass ceramic.

The high translucency and warm shade reproduction of VITA VM 11 combined with the opalescent effect of lithium silicate-type material facilitate highly esthetic results with a vivid play of colors.

Please note

* Can also be used for CELTRA® and CELTRA® CAD FC. CELTRA is a registered trademark of DeguDent GmbH.

Features/Benefits

- Natural shade and light effects thanks to the fine structure
- Excellent bonding ensures reliable results
- Easy handling for fast, highly-esthetic results
- Outstanding firing properties with dimensional stability retained even after repeated firing
- Excellent modeling characteristics for fast and accurate application of the ceramic
- Efficient and cost-effective processing
- Can be individualized using VITA AKZENT Plus and VITA INTERNO
- Process reliability thanks to long-standing experience
- Available in VITA SYSTEM 3D-Master and VITA classical A1 – D4 shades

VITAVM®11 kits

VITAVM®11 VITA classical A1-D4®



VITAVM®11 DENTINE KIT VITA classical A1-D4® BV11DKC

Quantity	Content	Material
8	12 g	VITA VM 11 TRANSPA DENTINE 0M1, A1 - A3.5, B2, C2, D2
2	12 g	VITA VM 11 ENAMEL ENL, END
1	50 ml	VITA LOW FUSING MODELLING LIQUID

VITAVM®11 VITA SYSTEM 3D-MASTER®



VITAVM®11 DENTINE KIT VITA SYSTEM 3D-MASTER® BV11DK

Quantity	Content	Material
6	12 g	VITA VM 11 TRANSPA DENTINE 0M1, 1M1 - 1M2, 2M2, 3M2, 4M2
2	12 g	VITA VM 11 ENAMEL ENL, END
1	50 ml	VITA LOW FUSING MODELLING LIQUID

VITAVM®11 additional powder

compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades



VITAVM®11 CREATIVE KIT BV11CKV1

Quantity	Content	Material
3	12 g	VITA VM 11 SUN DENTINE SD1 - SD3
2	12 g	VITA VM 11 MAMELON MM1, MM3
7	12 g	VITA VM 11 EFFECT ENAMEL EE1, EE3, EE5, EE7 - EE9, EE11
3	12 g	VITA VM 11 EFFECT CHROMA EC1, EC5, EC11
1	12 g	VITA VM 11 EFFECT PEARL EP1
4	12 g	VITA VM 11 EFFECT OPAL EO1 - EO3, EO5
2	12 g	VITA VM 11 ENAMEL ENL, END
1	12 g	VITA VM 11 NEUTRAL NT
1	12 g	VITA VM 11 WINDOW WIN
1	50 ml	VITA LOW FUSING MODELLING LIQUID
1		Shade guide for VITA VM 11 No. 1
1		Shade guide for VITA VM 11 No. 2

VITAVM®11 products

VITAVM®11 VITA classical A1-D4®

VITA VM 11 TRANSPA DENTINE, 12 g

Description	Powder no.	Prod. no. 12 g
A1	ENL	B4945112
A2	ENL	B4945212
A3	ENL	B4945312
A3.5	END	B4945412
B2	END	B4945712
C2	END	B4946112
D2	END	B4946412

VITAVM®11 VITA SYSTEM 3D-MASTER®

VITA VM 11 TRANSPA DENTINE, 12 g

Description	Powder no.	Prod. no. 12 g
0M1	ENL	B4906112
1M1	ENL	B4906412
1M2	ENL	B4906512
2M2	ENL	B4906912
3M2	ENL	B4907612
4M2	END	B4908312

VITAVM®11 additional powder

compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades

VITA VM 11 EFFECT ENAMEL, 12 g

Description	Prod. no. 12 g
EE1	B4909112
EE3	B4909312
EE5	B4909512
EE7	B4909712
EE8	B4909812
EE9	B4909912
EE11	B4910112

VITA VM 11 EFFECT CHROMA, 12 g

Description	Prod. no. 12 g
EC1	B4912112
EC5	B4912512
EC11	B4913112

VITA VM 11 EFFECT PEARL, 12 g

Description	Prod. no. 12 g
EP1	B4915112

VITA VM 11 MAMELON, 12 g

Description	Prod. no. 12 g
MM1	B4916112
MM3	B4916312

VITA VM 11 EFFECT OPAL, 12 g

Description	Prod. no. 12 g
E01	B4917112
E02	B4917212
E03	B4917312
E05	B4917512

VITA VM 11 WINDOW, 12 g

Description	Prod. no. 12 g
WIN	B4918112

VITA VM 11 ENAMEL, 12 g

Description	Prod. no. 12 g
ENL	B4919112
END	B4919212

VITA VM 11 NEUTRAL, 12 g

Description	Prod. no. 12 g
NT	B4920112

VITA VM 11 SUN DENTINE, 12 g

Description	Prod. no. 12 g
SD1	B4936112
SD2	B4936212
SD3	B4936312

VITAVM®9

Natural fine-structure feldspar veneering ceramic for zirconia frameworks, VITABLOCS Mark II and PM 9



Product description

VITA VM 9 is a highly-esthetic, fine-structure feldspar ceramic that is perfectly adapted to the CTE value of zirconia frameworks (approx. 10.5, e.g. such as VITA YZ) – regardless of whether they are white, colored or translucent.

In addition, you can also use VITA VM 9 for individualizing restorations made from VITABLOCS MARK II and VITA PM 9 press ceramic pellets.

The main components are pure-grade potash and albite feldspar materials that offer brilliant shade effects as well as optimum physical properties, such as extreme flexural strength values.

Thanks to its homogeneous, compact surface, VITA VM 9 supports outstanding grinding and polishing, particularly in situ. This results in smooth and densely compact surfaces.

Accumulation of plaque on the ceramic surface is significantly reduced, supporting the patient in caring for their high-quality restoration.

VITA VM 9 is a high-quality system component provided by the innovative VITA VM concept – the universal solution for all common indications. Metal-ceramics, all-ceramics or composite: the possibilities are endless!

VITA AKZENT Plus and VITA INTERNO stains are also available for individual finalization.

Features/Benefits

- Natural shade effects, light dynamics and optimum physical properties thanks to the fine structure
- Excellent bonding for reliable results and long-term success
- Minimum shrinkage behavior for accurate firing
- Excellent modeling characteristics for fast and accurate application of the ceramic
- Comprehensive range of additional materials for excellent photo-optic effects
- Can be individualized using VITA AKZENT Plus and VITA INTERNO
- Process reliability thanks to long-standing experience
- Available in VITA SYSTEM 3D-Master and VITA classical A1 – D4 shades

VITAVM®9 kits

VITAVM®9 VITA classical A1-D4®



VITAVM®9 BASIC KIT VITA classical A1-D4®

Basic kit for BASIC layering

BV9BK12CV2

Quantity	Content	Material
16	12 g	VITA VM 9 BASE DENTINE A1 - D4
5	12 g	VITA VM 9 CHROMA PLUS CP1 - CP5
2	12 g	VITA VM 9 ENAMEL ENL, END
1	12 g	VITA VM 9 NEUTRAL NT
1	12 g	VITA VM 9 WINDOW WIN
3	12 g	VITA VM 9 CORRECTIVE COR1 - COR3
3	12 g	VITA VM 9 SUN DENTINE SD1 - SD3
1	50 ml	VITA MODELLING FLUID RS
1		VITA Modisol®
1	2 pcs.	Fibrous pad firing support
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 9 classical
1		VITA classical A1-D4 shade guide



VITAVM®9 BASIC KIT SMALL VITA classical A1-D4®

Small basic kit for BASIC layering

BV9BKS12CV2

Quantity	Content	Material
6	12 g	VITA VM 9 BASE DENTINE A1 - A3.5, B3, D3
5	12 g	VITA VM 9 CHROMA PLUS CP1 - CP5
2	12 g	VITA VM 9 ENAMEL ENL, END
1	12 g	VITA VM 9 NEUTRAL NT
1	12 g	VITA VM 9 WINDOW WIN
1	12 g	VITA VM 9 CORRECTIVE COR1
3	12 g	VITA VM 9 SUN DENTINE SD1 - SD3
1	50 ml	VITA MODELLING FLUID RS
1		VITA Modisol®
1	2 pcs.	Fibrous pad firing support
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 9 classical
1		VITA classical A1-D4 shade guide



VITAVM®9 BUILD UP KIT VITA classical A1-D4®
Build-up kit for BUILD UP layering
BV9BUK12CV1

Quantity	Content	Material
16	12 g	VITA VM 9 TRANSPA DENTINE A1 - D4
1	50 ml	VITA MODELLING FLUID RS



VITAVM®9 BUILD UP KIT SMALL VITA classical A1-D4®
Small build-up kit for BUILD UP layering
BV9BUKS12CV1

Quantity	Content	Material
6	12 g	VITA VM 9 TRANSPA DENTINE A1 - A3.5, B3, D3
1	50 ml	VITA MODELLING FLUID RS



VITAVM®9 CLASSICAL COLOR KIT
Add-on kit for VITA VM 9 3D-MASTER users
BV9CCK12V2*

Quantity	Content	Material
16	12 g	VITA VM 9 BASE DENTINE A1 - D4
16	12 g	VITA VM 9 TRANSPA DENTINE A1 - D4
2	12 g	VITA VM 9 CHROMA PLUS CP1, CP5
1	50 ml	VITA MODELLING FLUID RS
1		Shade indicator for VITA VM 9 classical
1		VITA classical A1-D4 shade guide

* Assortment for VITA VM 9 SYSTEM 3D-MASTER customers who would like to add VITA classical A1-D4 to their existing range of materials



VITAVM®9 ONE COLOR KIT A2
One-shade kit for initial tests
BV9OCKA2V4

Quantity	Content	Material
1	12 g	VITA VM 9 CHROMA PLUS CP2
1	12 g	VITA VM 9 BASE DENTINE A2
1	12 g	VITA VM 9 ENAMEL ENL
1	12 g	VITA VM 9 WINDOW WIN
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	20 ml	VITA AKZENT Plus POWDER FLUID

VITAVM®9 VITA SYSTEM 3D-MASTER®



VITAVM®9 BASIC KIT VITA SYSTEM 3D-MASTER®

Basic kit for BASIC layering

BV9BK12V3

Quantity	Content	Material
26	12 g	VITA VM 9 BASE DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
2	12 g	VITA VM 9 ENAMEL ENL, END
1	12 g	VITA VM 9 NEUTRAL NT
1	12 g	VITA VM 9 WINDOW WIN
3	12 g	VITA VM 9 CORRECTIVE COR1 - COR3
3	12 g	VITA VM 9 CHROMA PLUS CP2 - CP4
3	12 g	VITA VM 9 SUN DENTINE SD1 - SD3
1	50 ml	VITA MODELLING FLUID RS
1		VITA Modisol®
1	2 pcs.	Fibrous pad firing support
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 9 3D-MASTER
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITAVM®9 BASIC KIT SMALL VITA SYSTEM 3D-MASTER®

Small basic kit for BASIC layering

BV9BKS12V3

Quantity	Content	Material
15	12 g	VITA VM 9 BASE DENTINE 1M1 - 4M3
2	12 g	VITA VM 9 ENAMEL ENL, END
1	12 g	VITA VM 9 NEUTRAL NT
1	12 g	VITA VM 9 WINDOW WIN
1	12 g	VITA VM 9 CORRECTIVE COR1
3	12 g	VITA VM 9 CHROMA PLUS CP2 - CP4
3	12 g	VITA VM 9 SUN DENTINE SD1 - SD3
1	50 ml	VITA MODELLING FLUID RS
1		VITA Modisol®
1	2 pcs.	Fibrous pad firing support
1		Firing tray G (pack of 2 trays and 20 pins)
1		Shade indicator for VITA VM 9 3D-MASTER
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITAVM®9 BLEACHED COLOR KIT VITA SYSTEM 3D-MASTER®
Very light shades to reproduce bleached teeth
BV9BLKV2

Quantity	Content	Material
3	12 g	VITA VM 9 BASE DENTINE 0M1 - 0M3
3	12 g	VITA VM 9 TRANSPA DENTINE 0M1 - 0M3
1	12 g	VITA VM 9 ENAMEL ENL
1	12 g	VITA VM 9 NEUTRAL NT
1	12 g	VITA VM 9 WINDOW WIN
1	50 ml	VITA MODELLING FLUID RS
1		VITA BLEACHED SHADE GUIDE, Shade Group 0



VITAVM®9 BUILD UP KIT VITA SYSTEM 3D-MASTER®
Build-up kit for BUILD UP layering
BV9BUK12V1

Quantity	Content	Material
26	12 g	VITA VM 9 TRANSPA DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
1	50 ml	VITA MODELLING FLUID RS

VITAVM®9 BUILD UP KIT SMALL VITA SYSTEM 3D-MASTER®
Small build-up kit for BUILD UP layering
BV9BUKS12V1

Quantity	Content	Material
15	12 g	VITA VM 9 TRANSPA DENTINE 1M1 - 4M3
1	50 ml	VITA MODELLING FLUID RS



VITAVM®9 ONE COLOR KIT 2M2
One-shade kit for initial tests
BV9OCK2M2V5

Quantity	Content	Material
1	12 g	VITA VM 9 CHROMA PLUS CP2
1	12 g	VITA VM 9 BASE DENTINE 2M2
1	12 g	VITA VM 9 ENAMEL ENL
1	12 g	VITA VM 9 WINDOW WIN
1	5 g	VITA AKZENT Plus GLAZE POWDER
1	50 ml	VITA MODELLING FLUID RS
1	20 ml	VITA AKZENT Plus POWDER FLUID

VITAVM®9 additional powder
compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades

VITAVM®9 ADD-ON KIT

Set including individualization materials and accessories
BV9AOKV2

Quantity	Content	Material
8	12 g	VITA VM 9 ADD-ON ADD1 - ADD8
1	5 g	VITA AKZENT Plus GLAZE LT POWDER
1	50 ml	VITA VM MODELLING LIQUID
1	20 ml	VITA AKZENT Plus POWDER FLUID
1	2 pcs.	Fibrous pad firing support
1		Shade guide for VITA VM 9 ADD-ON
1		Paint brush no. 3/0 synthetic



VITAVM®9 ESTHETIC KIT for VITABLOCS for CEREC

Additional materials for customizing VITABLOCS
BV9EKC2

Quantity	Content	Material
1		VITABLOCS Mark II Sample Set VITA SYSTEM 3D-MASTER
1	12 g	VITA VM 9 WINDOW WIN
1	12 g	VITA VM 9 NEUTRAL NT
2	12 g	VITA VM 9 ENAMEL ENL, END
1	12 g	VITA VM 9 EFFECT PEARL EP1
2	12 g	VITA VM 9 EFFECT ENAMEL EE1, EE10
1	12 g	VITA VM 9 CORRECTIVE COR1
1	4 g	VITA AKZENT Plus FINISHING AGENT PASTE
1	4 g	VITA AKZENT Plus GLAZE PASTE
1	12 g	VITA VM 9 EFFECT OPAL EO2
2	12 g	VITA VM 9 EFFECT CHROMA EC1, EC4
1	12 g	VITA VM 9 MAMELON MM2
1	50 ml	VITA MODELLING FLUID RS
1	15 ml	VITA AKZENT Plus PASTE FLUID
1		Firing tray W for crowns and bridges
1		2 Platinum pins (1 anterior and 1 posterior pin)



VITAVM®9 GINGIVA KIT
Natural-looking gingiva pastes
BV9GKV1

Quantity	Content	Material
5	12 g	VITA VM 9 GINGIVA G1 - G5
1		Shade guide for VITA VM GINGIVA / MARGIN



VITAVM®9 MARGIN KIT
Only for minor corrections in the area of margins
BV9MKV1

Quantity	Content	Material
6	12 g	VITA VM 9 MARGIN M1, M4 - M5, M7 - M9
1		Shade guide for VITA VM GINGIVA / MARGIN



VITAVM®9 PROFESSIONAL KIT SMALL
Additional materials for achieving special shade effects and individualizing the restoration
BV9PKS

Quantity	Content	Material
5	12 g	VITA VM 9 EFFECT CHROMA EC1, EC4, EC6, EC8 - EC9
1	12 g	VITA VM 9 MAMELON MM2
1	12 g	VITA VM 9 EFFECT PEARL EP1
1	12 g	VITA VM 9 EFFECT OPAL EO2
7	12 g	VITA VM 9 EFFECT ENAMEL EE1, EE3, EE7 - EE11
1		Shade guide for VITA VM PROFESSIONAL KIT SMALL



VITAVM®9 PROFESSIONAL KIT
Additional materials for achieving special shade effects and individualizing the restoration
BV9PKV2

Quantity	Content	Material
11	12 g	VITA VM 9 EFFECT CHROMA EC1 - EC11
3	12 g	VITA VM 9 MAMELON MM1 - MM3
3	12 g	VITA VM 9 EFFECT PEARL EP1 - EP3
5	12 g	VITA VM 9 EFFECT OPAL EO1 - EO5
11	12 g	VITA VM 9 EFFECT ENAMEL EE1 - EE11
6	12 g	VITA VM 9 EFFECT LINER EL1 - EL6
1		Shade guide for VITA VM EFFECT CHROMA
1		Shade guide for VITA VM EFFECT LINER / EFFECT OPAL / MAMELON
1		Shade guide for VITA VM EFFECT ENAMEL / EFFECT PEARL

VITAVM®9 products**VITAVM®9 VITA classical A1-D4®**

VITA VM 9 BASE DENTINE, 12 g / 50 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
A1	ENL	B4243112	B4243150
A2	ENL	B4243212	B4243250
A3	ENL	B4243312	B4243350
A3.5	END	B4243412	B4243450
A4	END	B4243512	B4243550
B1	END	B4243612	B4243650
B2	END	B4243712	B4243750
B3	END	B4243812	B4243850
B4	END	B4243912	B4243950
C1	END	B4244012	B4244050
C2	END	B4244112	B4244150
C3	ENL	B4244212	B4244250
C4	ENL	B4244312	B4244350
D2	END	B4244412	B4244450
D3	END	B4244512	B4244550
D4	END	B4244612	B4244650

VITA VM 9 TRANSPA DENTINE, 12 g / 50 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
A1	ENL	B4245112	B4245150
A2	ENL	B4245212	B4245250
A3	ENL	B4245312	B4245350
A3.5	END	B4245412	B4245450
A4	END	B4245512	B4245550
B1	END	B4245612	B4245650
B2	END	B4245712	B4245750
B3	END	B4245812	B4245850
B4	END	B4245912	B4245950
C1	END	B4246012	B4246050
C2	END	B4246112	B4246150
C3	ENL	B4246212	B4246250
C4	ENL	B4246312	B4246350
D2	END	B4246412	B4246450
D3	END	B4246512	B4246550
D4	END	B4246612	B4246650

VITAVM® 9 VITA SYSTEM 3D-MASTER®

VITA VM 9 BASE DENTINE, 12 g / 50 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
0M1	ENL	B4203112	B4203150
0M2	ENL	B4203212	B4203250
0M3	ENL	B4203312	B4203350
1M1	ENL	B4203412	B4203450
1M2	ENL	B4203512	B4203550
2L1.5	ENL	B4203612	B4203650
2L2.5	ENL	B4203712	B4203750
2M1	ENL	B4203812	B4203850
2M2	ENL	B4203912	B4203950
2M3	ENL	B4204012	B4204050
2R1.5	ENL	B4204112	B4204150
2R2.5	ENL	B4204212	B4204250
3L1.5	ENL	B4204312	B4204350
3L2.5	ENL	B4204412	B4204450
3M1	ENL	B4204512	B4204550
3M2	ENL	B4204612	B4204650
3M3	ENL	B4204712	B4204750

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
3R1.5	ENL	B4204812	B4204850
3R2.5	ENL	B4204912	B4204950
4L1.5	END	B4205012	B4205050
4L2.5	END	B4205112	B4205150
4M1	END	B4205212	B4205250
4M2	END	B4205312	B4205350
4M3	END	B4205412	B4205450
4R1.5	END	B4205512	B4205550
4R2.5	END	B4205612	B4205650
5M1	END	B4205712	B4205750
5M2	END	B4205812	B4205850
5M3	END	B4205912	B4205950

VITA VM 9 TRANSPA DENTINE, 12 g / 50 g

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
0M1	ENL	B4206112	B4206150
0M2	ENL	B4206212	B4206250
0M3	ENL	B4206312	B4206350
1M1	ENL	B4206412	B4206450
1M2	ENL	B4206512	B4206550
2L1.5	ENL	B4206612	B4206650
2L2.5	ENL	B4206712	B4206750
2M1	ENL	B4206812	B4206850
2M2	ENL	B4206912	B4206950
2M3	ENL	B4207012	B4207050
2R1.5	ENL	B4207112	B4207150
2R2.5	ENL	B4207212	B4207250
3L1.5	ENL	B4207312	B4207350
3L2.5	ENL	B4207412	B4207450
3M1	ENL	B4207512	B4207550
3M2	ENL	B4207612	B4207650
3M3	ENL	B4207712	B4207750
3R1.5	ENL	B4207812	B4207850
3R2.5	ENL	B4207912	B4207950
4L1.5	END	B4208012	B4208050
4L2.5	END	B4208112	B4208150
4M1	END	B4208212	B4208250
4M2	END	B4208312	B4208350
4M3	END	B4208412	B4208450

Description	Powder no.	Prod. no. 12 g	Prod. no. 50 g
4R1.5	END	B4208512	B4208550
4R2.5	END	B4208612	B4208650
5M1	END	B4208712	B4208750
5M2	END	B4208812	B4208850
5M3	END	B4208912	B4208950

VITAVM®9 additional powder
compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades

VITA VM 9 EFFECT ENAMEL, 12 g

Description	Prod. no. 12 g
EE1	B4209112
EE2	B4209212
EE3	B4209312
EE4	B4209412
EE5	B4209512
EE6	B4209612
EE7	B4209712
EE8	B4209812
EE9	B4209912
EE10	B4210012
EE11	B4210112

VITA VM 9 EFFECT CHROMA, 12 g

Description	Prod. no. 12 g
EC1	B4212112
EC2	B4212212
EC3	B4212312
EC4	B4212412
EC5	B4212512
EC6	B4212612
EC7	B4212712
EC8	B4212812
EC9	B4212912
EC10	B4213012
EC11	B4213112

VITA VM 9 EFFECT LINER, 12 g

Description	Prod. no. 12 g
EL1	B4214112
EL2	B4214212
EL3	B4214312
EL4	B4214412
EL5	B4214512
EL6	B4214612

VITA VM 9 EFFECT PEARL, 12 g

Description	Prod. no. 12 g
EP1	B4215112
EP2	B4215212
EP3	B4215312

VITA VM 9 MAMELON, 12 g

Description	Prod. no. 12 g
MM1	B4216112
MM2	B4216212
MM3	B4216312

VITA VM 9 EFFECT OPAL, 12 g

Description	Prod. no. 12 g
E01	B4217112
E02	B4217212
E03	B4217312
E04	B4217412
E05	B4217512

VITA VM 9 WINDOW, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
WIN	B4218112	B4218150	B42181250

VITA VM 9 ENAMEL, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
ENL	B4219112	B4219150	B42191250
END	B4219212	B4219250	B42192250

VITA VM 9 NEUTRAL, 12 g / 50 g / 250 g

Description	Prod. no. 12 g	Prod. no. 50 g	Prod. no. 250 g
NT	B4220112	B4220150	B42201250

VITA VM 9 EFFECT BONDER, 12 g

Description	Prod. no. 12 g
EB0	B4221012
EB1	B4221112
EB2	B4221212
EB3	B4221312
EB4	B4221412
EB5	B4221512

VITA VM 9 CORRECTIVE, 12 g

Description	Prod. no. 12 g
COR1	B4222112
COR2	B4222212
COR3	B4222312

VITA VM 9 GINGIVA, 12 g

Description	Prod. no. 12 g
G1	B4223112
G2	B4223212
G3	B4223312
G4	B4223412
G5	B4223512

VITA VM 9 CHROMA PLUS, 12 g

Description	Prod. no. 12 g
CP1	B4224112
CP2	B4224212
CP3	B4224312
CP4	B4224412
CP5	B4224512

VITA VM 9 MARGIN, 12 g

Description	Prod. no. 12 g
M1	B4226112
M4	B4226412
M5	B4226512
M7	B4226712
M8	B4226812
M9	B4226912

VITA VM 9 ADD-ON, 12 g

Description	Prod. no. 12 g
ADD1	B4229112
ADD2	B4229212
ADD3	B4229312
ADD4	B4229412
ADD5	B4229512
ADD6	B4229612
ADD7	B4229712
ADD8	B4229812

VITA VM 9 SUN DENTINE, 12 g

Description	Prod. no. 12 g
SD1	B4236112
SD2	B4236212
SD3	B4236312

VITAVM® LC

Light-curing, microparticle composite for extraoral use with fixed and removable restorations



Product description

VITAVM® LC & VITAVM® LC flow

Simple to use. Cost-effective and esthetic results. Versatile range of applications: the VITA VM LC product family meets all your expectations when it comes to a veneering composite for the extraoral fabrication of restorations.

Whether it's in the form of a paste composite, the low-viscosity flow materials or freely combined, VITA VM LC offers brilliant results. VITA VM LC is also incredibly similar to ceramic.

Features/Benefits

- Two forms of processing for almost unlimited possibilities
- Perfectly tailored to suit any indication
- Brilliant esthetics with excellent shade reproduction based on the original VITA shade formula
- Outstanding material properties
- System-based, modular product range for efficient design options

Indications

Indication:

- Full and partial veneering of metal frameworks: crowns, bridges, telescopic crowns, implant suprastructures
- Inlays, veneers

Application areas:

- Individualization and layering-over VITA ENAMIC
- Veneering partially yttrium-stabilized ZrO₂ frameworks (CTE approx. 10.5 · 10⁻⁶ · K⁻¹), such as VITA YZ
- Individualization of VITA denture teeth
- Reproduction of gingival areas
- Veneering of removable and partially removable dentures (according to the manufacturer's information) made of
 - polyether ether ketone (PEEK) with a filler content of up to 20%, such as BioHPP/Bredent
 - PEEK-OPTIMA® LT1 polymer, such as, Juvora, InnoBlanc Medical

Long-term temporaries:

- Individualizing and layering-over long-term temporaries made from VITA CAD-Temp
- Metal-free crowns and three-unit anterior bridges made from VITA VM LC

Shades

Available in VITA SYSTEM 3D-MASTER® and VITA classical A1–D4® shades

VITAVM® LC kits

VITAVM® LC VITA classical A1-D4®



VITAVM®LC combi BASIC KIT VITA classical A1-D4®

Basic kit for BASIC layering (DENTINE paste / ENAMEL & additional powder flow)

CVLCCBKC

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
15	3 g	VITA VM LC OPAQUE PASTE A1 - D4
15	4 g	VITA VM LC BASE DENTINE A1 - D4
2	4 g	VITA VM LC flow ENAMEL ENL, END
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
1	20 ml	VITA VM LC GEL
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA classical A1-D4 shade guide



VITAVM®LC combi BASIC KIT SMALL VITA classical A1-D4®
Small basic kit for BASIC layering (DENTINE paste / ENAMEL & additional powder flow)
CVLCCBKSC

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
6	3 g	VITA VM LC OPAQUE PASTE A1 - A3.5, B3, D3
6	4 g	VITA VM LC BASE DENTINE A1 - A3.5, B3, D3
4	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP3, CP5
1	20 ml	VITA VM LC GEL
2	4 g	VITA VM LC flow ENAMEL ENL, END
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA classical A1-D4 shade guide



VITAVM®LC combi TEST KIT A3
One-shade kit for initial tests (DENTINE paste / ENAMEL & additional powder flow)
CVLCCTKA3

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
1	3 g	VITA VM LC OPAQUE PASTE A3
1	4 g	VITA VM LC BASE DENTINE A3
1	4 g	VITA VM LC flow ENAMEL ENL
1	3 g	VITA VM LC flow CHROMA PLUS CP2
1	5 ml	VITA VM LC Gel
1	10 ml	VITA VM LC MODELLING LIQUID
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. 1 correction brush, red marten
1		VITA shade indicator stick A3



VITAVM®LC flow BASIC KIT VITA classical A1-D4®

Basic kit for BASIC layering DENTINE, ENAMEL additional powder flow
CVLCFBKCV2

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
15	3 g	VITA VM LC OPAQUE PASTE A1 - D4
15	4 g	VITA VM LC flow BASE DENTINE A1 - D4
2	4 g	VITA VM LC flow ENAMEL ENL, END
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
1	50 ml	VITA VM LC CLEANER
1	20 ml	VITA VM LC GEL
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA classical A1-D4 shade guide

VITAVM®LC flow BASIC KIT SMALL VITA classical A1-D4®
Small basic kit for BASIC layering DENTINE, ENAMEL additional powder flow
CVLCFBKSCV2

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
6	3 g	VITA VM LC OPAQUE PASTE A1 - A3.5, B3, D3
6	4 g	VITA VM LC flow BASE DENTINE A1 - A3.5, B3, D3
2	4 g	VITA VM LC flow ENAMEL ENL, END
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
4	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP3, CP5
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1	20 ml	VITA VM LC GEL
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA classical A1-D4 shade guide



VITAVM®LC flow TEST KIT A3
One-shade kit for initial tests (DENTINE, ENAMEL & additional powder flow)
CVLCFTKA3

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
1	3 g	VITA VM LC OPAQUE PASTE A3
1	4 g	VITA VM LC flow BASE DENTINE A3
1	4 g	VITA VM LC flow ENAMEL ENL
1	5 ml	VITA VM LC Gel
1	10 ml	VITA VM LC MODELLING LIQUID
1	3 g	VITA VM LC flow CHROMA PLUS CP2
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. 1 correction brush, red marten
1		VITA shade indicator stick A3



VITAVM®LC paste BASIC KIT VITA classical A1-D4®

Basic kit for BASIC layering (DENTINE, ENAMEL additional powder paste)
CVLCPBKC

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
15	3 g	VITA VM LC OPAQUE PASTE A1 - D4
15	4 g	VITA VM LC BASE DENTINE A1 - D4
2	4 g	VITA VM LC ENAMEL ENL, END
1	4 g	VITA VM LC NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
1	20 ml	VITA VM LC GEL
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA classical A1-D4 shade guide



VITAVM®LC paste BASIC KIT SMALL VITA classical A1-D4®
Small basic kit for BASIC layering DENTINE, ENAMEL additional powder
paste
CVLCPBKSC

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
6	3 g	VITA VM LC OPAQUE PASTE A1 - A3.5, B3, D3
6	4 g	VITA VM LC BASE DENTINE A1 - A3.5, B3, D3
4	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP3, CP5
1	20 ml	VITA VM LC GEL
2	4 g	VITA VM LC ENAMEL ENL, END
1	4 g	VITA VM LC NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA classical A1-D4 shade guide



VITAVM®LC paste TEST KIT A3
One-shade kit for initial tests (DENTINE, ENAMEL & additional powder
paste)
CVLCPTKA3

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
1	3 g	VITA VM LC OPAQUE PASTE A3
1	4 g	VITA VM LC BASE DENTINE A3
1	4 g	VITA VM LC ENAMEL ENL
1	4 g	VITA VM LC CHROMA PLUS CP2
1	5 ml	VITA VM LC Gel
1	10 ml	VITA VM LC MODELLING LIQUID
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		VITA shade indicator stick A3

VITAVM® LC VITA SYSTEM 3D-MASTER®



VITAVM® LC combi BASIC KIT VITA SYSTEM 3D-MASTER®

Basic kit for BASIC layering (DENTINE paste / ENAMEL & additional powder flow)

CVLCCBK

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
15	10 g	VITA VM LC OPAQUE 2L2.5 - 5M3
11	3 g	VITA VM LC OPAQUE PASTE 1M1 - 1M2, 2L1.5, 2M1 - 2M3, 3L1.5, 3M2 - 3M3, 3R2.5, 4M2
26	4 g	VITA VM LC BASE DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
2	4 g	VITA VM LC flow ENAMEL ENL, END
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
2	5 ml	VITA VM LC OPAQUE LIQUID
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1	20 ml	VITA VM LC GEL
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		Measuring spoon, 0,1 ml / Mixing spatula 2-sided
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITAVM®LC combi BASIC KIT SMALL VITA SYSTEM 3D-MASTER®
Small basic kit for BASIC layering (DENTINE paste / ENAMEL & additional powder flow)
CVLCCBKS

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
11	3 g	VITA VM LC OPAQUE PASTE 1M1 - 1M2, 2L1.5, 2M1 - 2M3, 3L1.5, 3M2 - 3M3, 3R2.5, 4M2
11	4 g	VITA VM LC BASE DENTINE 1M1 - 1M2, 2L1.5, 2M1 - 2M3, 3L1.5, 3M2 - 3M3, 3R2.5, 4M2
2	4 g	VITA VM LC flow ENAMEL ENL, END
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1	20 ml	VITA VM LC GEL
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITAVM®LC combi TEST KIT 2M2
One-shade kit for initial tests (DENTINE paste / ENAMEL & additional powder flow)
CVLCCTK2M2

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
1	3 g	VITA VM LC OPAQUE PASTE 2M2
1	4 g	VITA VM LC BASE DENTINE 2M2
1	4 g	VITA VM LC flow ENAMEL ENL
1	3 g	VITA VM LC flow CHROMA PLUS CP2
1	5 ml	VITA VM LC Gel
1	10 ml	VITA VM LC MODELLING LIQUID
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		VITA shade indicator stick 2M2
1		Brush No. 1 correction brush, red marten

VITAVM®LC flow BASIC KIT SMALL VITA SYSTEM 3D-MASTER®
Small basic kit for BASIC layering (DENTINE, ENAMEL additional powder flow)
CVLCFBKSV1

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
8	3 g	VITA VM LC OPAQUE PASTE 0M1, 1M1, 2M1 - 2M3, 3M2 - 3M3, 3R2.5
8	4 g	VITA VM LC flow BASE DENTINE 0M1, 1M1, 2M1 - 2M3, 3M2 - 3M3, 3R2.5
1	4 g	VITA VM LC flow ENAMEL ENL
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1	20 ml	VITA VM LC GEL
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA Toothguide 3D-Master with Bleached Shade Guide

VITAVM®LC flow TEST KIT 2M2
One-shade kit for initial tests (DENTINE, ENAMEL & additional powder flow)
CVLCFTK2M2

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
1	3 g	VITA VM LC OPAQUE PASTE 2M2
1	4 g	VITA VM LC flow BASE DENTINE 2M2
1	4 g	VITA VM LC flow ENAMEL ENL
1	3 g	VITA VM LC flow CHROMA PLUS CP2
1	5 ml	VITA VM LC Gel
1	10 ml	VITA VM LC MODELLING LIQUID
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. 1 correction brush, red marten
1		VITA shade indicator stick 2M2



VITAVM®LC paste BASIC KIT VITA SYSTEM 3D-MASTER®
Basic kit for BASIC layering DENTINE, ENAMEL additional powder paste
CVLCPBK

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
15	10 g	VITA VM LC OPAQUE 2L2.5 - 5M3
11	3 g	VITA VM LC OPAQUE PASTE 1M1 - 1M2, 2L1.5, 2M1 - 2M3, 3L1.5, 3M2 - 3M3, 3R2.5, 4M2
26	4 g	VITA VM LC BASE DENTINE 1M1 - 1M2, 2L1.5 - 2L2.5, 2M1 - 2M3, 2R1.5 - 2R2.5, 3L1.5 - 3L2.5, 3M1 - 3M3, 3R1.5 - 3R2.5, 4L1.5 - 4L2.5, 4M1 - 4M3, 4R1.5 - 4R2.5, 5M1 - 5M3
2	4 g	VITA VM LC ENAMEL ENL, END
1	4 g	VITA VM LC NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
2	5 ml	VITA VM LC OPAQUE LIQUID
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1	20 ml	VITA VM LC GEL
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Mixing spatula, Acrylic for PASTE OPAQUE
1		Measuring spoon, 0,1 ml / Mixing spatula 2-sided
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITAVM®LC paste BASIC KIT SMALL VITA SYSTEM 3D-MASTER®
Small basic kit for BASIC layering DENTINE, ENAMEL additional powder
paste
CVLCPBKS

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
11	3 g	VITA VM LC OPAQUE PASTE 1M1 - 1M2, 2L1.5, 2M1 - 2M3, 3L1.5, 3M2 - 3M3, 3R2.5, 4M2
11	4 g	VITA VM LC BASE DENTINE 1M1 - 1M2, 2L1.5, 2M1 - 2M3, 3L1.5, 3M2 - 3M3, 3R2.5, 4M2
2	4 g	VITA VM LC ENAMEL ENL, END
1	4 g	VITA VM LC NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
1	10 ml	VITA VM LC MODELLING LIQUID
1	30 ml	VITA VM LC SEPARATOR
1	50 ml	VITA VM LC CLEANER
1	20 ml	VITA VM LC GEL
1		VITA Mixing palette black porcelain with lid, 8,5 x 11 cm
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		Brush No. E 0 staining
1		Shade guide for VITA VM INDIVIDUAL
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITAVM®LC paste TEST KIT 2M2
One-shade kit for initial tests (DENTINE, ENAMEL & additional powder
paste)
CVLCPTK2M2

Quantity	Content	Material
1	3 g	VITA VM LC PRE OPAQUE PRE
1	3 g	VITA VM LC OPAQUE PASTE 2M2
1	4 g	VITA VM LC BASE DENTINE 2M2
1	4 g	VITA VM LC CHROMA PLUS CP2
1	5 ml	VITA VM LC Gel
1	10 ml	VITA VM LC MODELLING LIQUID
1		White brush-holder
1	50 pcs.	Disposable brush heads
1		VITA shade indicator stick 2M2
1	4 g	VITA VM LC ENAMEL ENL

VITAVM®LC additional powder
compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades



VITAVM®LC flow CHROMA KIT
Shade intensification kit for the neck and body area
CVLCFCHK

Quantity	Content	Material
5	3 g	VITA VM LC flow CHROMA PLUS CP1 - CP5
1		Shade guide for VITA VM INDIVIDUAL



VITAVM®LC flow CREATIVE KIT
for the layering of CAD-Temp and for individualizing VITA denture teeth
CVLCFCKV1

Quantity	Content	Material
4	2 g	VITA VM LC PAINT PT1, PT5, PT15, PT17
1	4 g	VITA VM LC flow ENAMEL ENL
2	3 g	VITA VM LC flow EFFECT ENAMEL EE6, EE9
1	4 g	VITA VM LC flow NEUTRAL NT
1	3 g	VITA VM LC flow WINDOW WIN
1	5 ml	VITA VM LC Gel
1	10 ml	VITA VM LC MODELLING LIQUID
1		Brush No. E 0 staining



VITAVM®LC flow GINGIVA KIT
Reconstruction kit for gingiva parts
CVLCFGK

Quantity	Content	Material
1	3 g	VITA VM LC GINGIVA OPAQUE PASTE GOL
5	3 g	VITA VM LC flow GINGIVA G1 - G5
1		Shade guide for VITA VM INDIVIDUAL



VITAVM®LC flow PROFESSIONAL KIT
Additional materials for achieving special shade effects and individualizing the restoration
CVLCFPKV1

Quantity	Content	Material
8	3 g	VITA VM LC flow EFFECT ENAMEL EE1 - EE2, EE5 - EE7, EE9, EE11 - EE12
1	3 g	VITA VM LC flow WINDOW WIN
1	4 g	VITA VM LC flow NEUTRAL NT
1		Shade guide for VITA VM INDIVIDUAL



VITAVM®LC PAINT KIT
Intensive shades
CVLCPAKV4

Quantity	Content	Material
10	2 g	VITA VM LC PAINT PT1, PT3, PT5, PT8 - PT9, PT12 - PT13, PT15, PT17, PT19
1	3 g	VITA VM LC flow WINDOW WIN
1		Shade guide for VITA VM INDIVIDUAL



VITAVM®LC paste GINGIVA KIT
Reconstruction kit for gingiva parts
CVLCPGK

Quantity	Content	Material
1	3 g	VITA VM LC GINGIVA OPAQUE PASTE GOL
4	4 g	VITA VM LC GINGIVA G1 - G2, G4 - G5
1		Shade guide for VITA VM INDIVIDUAL

VITAVM® LC products

VITAVM® LC VITA classical A1-D4®

VITA VM LC BASE DENTINE, 4 g

Description	Powder no.	Prod. no. 4 g
A1	ENL	C444314
A2	ENL	C444324
A3	ENL	C444334
A3.5	END	C444344
A4	END	C444354
B2	ENL	C444374
B3	END	C444384
B4	END	C444394
C1	END	C444404
C2	ENL	C444414
C3	END	C444424
C4	END	C444434
D2	END	C444444
D3	END	C444454
D4	END	C444464

VITA VM LC flow BASE DENTINE, 4 g

Description	Powder no.	Prod. no. 4 g
A1	ENL	C44431F4
A2	ENL	C44432F4
A3	ENL	C44433F4
A3.5	END	C44434F4
A4	END	C44435F4
B2	END	C44437F4
B3	END	C44438F4
B4	END	C44439F4
C1	END	C44440F4
C2	ENL	C44441F4
C3	END	C44442F4
C4	END	C44443F4
D2	END	C44444F4
D3	END	C44445F4
D4	END	C44446F4

VITA VM LC OPAQUE PASTE, 3 g

Description	Prod. no. 3 g
A1	C443013
A2	C443023
A3	C443033
A3.5	C443043
A4	C443053
B2	C443073
B3	C443083
B4	C443093
C1	C443103
C2	C443113
C3	C443123
C4	C443133
D2	C443143
D3	C443153
D4	C443163

VITA VM LC OPAQUE, 10g

Description	Prod. no. 10g
A1	C4440110
A2	C4440210
A3	C4440310
A3.5	C4440410
B3	C4440810
D3	C4441510

VITAVM® LC VITA SYSTEM 3D-MASTER®

VITA VM LC BASE DENTINE, 4 g

Description	Powder no.	Prod. no. 4 g
0M1	ENL	C440314
1M1	ENL	C440344
1M2	ENL	C440354
2L1.5	ENL	C440364
2L2.5	ENL	C440374
2M1	ENL	C440384
2M2	ENL	C440394
2M3	ENL	C440404
2R1.5	ENL	C440414

Description	Powder no.	Prod. no. 4 g
2R2.5	ENL	C440424
3L1.5	ENL	C440434
3L2.5	ENL	C440444
3M1	ENL	C440454
3M2	ENL	C440464
3M3	ENL	C440474
3R1.5	ENL	C440484
3R2.5	ENL	C440494
4L1.5	END	C440504
4L2.5	END	C440514
4M1	END	C440524
4M2	END	C440534
4M3	END	C440544
4R1.5	END	C440554
4R2.5	END	C440564
5M1	END	C440574
5M2	END	C440584
5M3	END	C440594

VITA VM LC flow BASE DENTINE, 4 g

Description	Powder no.	Prod. no. 4 g
0M1	ENL	C44031F4
1M1	ENL	C44034F4
2M1	ENL	C44038F4
2M2	ENL	C44039F4
2M3	ENL	C44040F4
3M2	ENL	C44046F4
3M3	ENL	C44047F4
3R2.5	ENL	C44049F4

VITA VM LC OPAQUE PASTE, 3 g

Description	Prod. no. 3 g
0M1	C445013
1M1	C445043
1M2	C445053
2L1.5	C445063
2M1	C445083
2M2	C445093
2M3	C445103
3L1.5	C445133
3M2	C445163
3M3	C445173
3R2.5	C445193
4M2	C445233

VITA VM LC OPAQUE, 10 g

Description	Prod. no. 10 g
0M1	C4400110
1M1	C4400410
1M2	C4400510
2L1.5	C4400610
2L2.5	C4400710
2M1	C4400810
2M2	C4400910
2M3	C4401010
2R1.5	C4401110
2R2.5	C4401210
3L1.5	C4401310
3L2.5	C4401410
3M1	C4401510
3M2	C4401610
3M3	C4401710
3R1.5	C4401810
3R2.5	C4401910
4L1.5	C4402010
4L2.5	C4402110
4M1	C4402210
4M2	C4402310
4M3	C4402410

Description	Prod. no. 10 g
4R1.5	C4402510
4R2.5	C4402610
5M1	C4402710
5M2	C4402810
5M3	C4402910

VITAVM® LC additional powder
compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades

VITA VM LC CHROMA PLUS, 4 g

Description	Prod. no. 4 g
CP1	C442414
CP2	C442424
CP3	C442434
CP4	C442444
CP5	C442454

VITA VM LC COLOR OPAQUE, 10 g

Description	Prod. no. 10 g
C01	C4428110
C02	C4428210
C03	C4428310

VITA VM LC EFFECT ENAMEL, 4 g

Description	Prod. no. 4 g
EE12	C441024
EE2	C440924
EE6	C440964
EE9	C440994

VITA VM LC ENAMEL, 4 g

Description	Prod. no. 4 g
END	C441924
ENL	C441914

VITA VM LC flow CHROMA PLUS, 3 g

Description	Prod. no. 3 g
CP1	C44241F3
CP2	C44242F3
CP3	C44243F3
CP4	C44244F3
CP5	C44245F3

VITA VM LC flow EFFECT ENAMEL, 3 g

Description	Prod. no. 3 g
EE1	C44091F3
EE11	C44101F3
EE12	C44102F3
EE2	C44092F3
EE5	C44095F3
EE6	C44096F3
EE7	C44097F3
EE9	C44099F3

VITA VM LC flow ENAMEL, 4 g

Description	Prod. no. 4 g
END	C44192F4
ENL	C44191F4

VITA VM LC flow GINGIVA, 3 g

Description	Prod. no. 3 g
G1	C44231F3
G2	C44232F3
G3	C44233F3
G4	C44234F3
G5	C44235F3

VITA VM LC flow NEUTRAL, 4 g

Description	Prod. no. 4 g
NT	C44201F4

VITA VM LC flow WINDOW, 3 g

Description	Prod. no. 3 g
WIN	C44181F3

VITA VM LC GINGIVA OPAQUE PASTE, 3 g

Description	Prod. no. 3 g
GOL	C443713

VITA VM LC GINGIVA OPAQUE, 10 g

Description	Prod. no. 10 g
GOL	C4427110

VITA VM LC GINGIVA, 4 g

Description	Prod. no. 4 g
G1	C442314
G2	C442324
G4	C442344
G5	C442354

VITA VM LC NEUTRAL, 4 g

Description	Prod. no. 4 g
NT	C442014

Description		Prod. no.
VITA VM LC OPAQUE LIQUID, 5 ml Opaque liquid		CVML0L5

VITA VM LC PAINT, 2 g

Description	Prod. no. 2 g
PT1	C444012
PT12	C444122
PT13	C444132
PT15	C444152
PT17	C444172
PT19	C444192
PT3	C444032
PT5	C444052
PT8	C444082
PT9	C444092

VITA VM LC PRE OPAQUE, 3 g

Description	Prod. no. 3 g
PRE	C445913

VITAVM[®]CC

Filler-free cold-curing polymer resin for extra-oral fabrication



Product description

VITA VM CC is a filler-free cold-curing resin. It is the ideal material for the fabrication of indirect temporary restorations.

VITA VM CC exhibits excellent processing characteristics. Processing is particularly simplified thanks to the exceptional mixing behavior and the superior flow properties.

Accumulation of plaque can be considerably reduced since the material can be easily polished and features a homogeneous structure. The smooth surface of the temporary restoration is pleasant to the patient and offers superior wearing comfort.

Features/Benefits

- Outstanding flowability
- First-rate mixing behavior
- Easy polishing
- Greater color stability than comparable systems
- The ideal combination of bending strength and elasticity modulus ensure
- extreme material stability and lower fracture rates
- Bubble-free curing, with or without pressure pot

Indications

- Fabrication of temporary crowns and bridges
- Fixation of VITA denture teeth, e.g. on CoCr structures
- Repairs
- Covering of attachments in combined restorations

Shades

Available in VITA SYSTEM 3D-MASTER and VITA classical A1–D4 shades

VITAVM® CC kits

VITAVM® CC VITA classical A1-D4®



VITAVM® CC BASIC KIT VITA classical A1-D4®

Basic kit for BASIC layering

CVCCBKC

Quantity	Content	Material
15	30 g	VITA VM CC BASE DENTINE A1 - D4
2	30 g	VITA VM CC ENAMEL ENL, END
1	30 g	VITA VM CC WINDOW WIN
1	100 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA classical A1-D4 shade guide



VITAVM® CC BASIC KIT SMALL VITA classical A1-D4®

Small basic kit for BASIC layering

CVCCBKSC

Quantity	Content	Material
6	30 g	VITA VM CC BASE DENTINE A1 - A3.5, B3, D3
2	30 g	VITA VM CC ENAMEL ENL, END
1	30 g	VITA VM CC WINDOW WIN
1	30 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA classical A1-D4 shade guide



VITA VM CC ONE COLOR KIT A2

One-shade kit for initial tests

CVCCOCKA2

Quantity	Content	Material
1	30 g	VITA VM CC BASE DENTINE A2
1	30 g	VITA VM CC ENAMEL ENL
1	30 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA shade indicator stick A2

VITA VM CC ONE COLOR KIT A3

One-shade kit for initial tests

CVCCOCKA3

Quantity	Content	Material
1	30 g	VITA VM CC BASE DENTINE A3
1	30 g	VITA VM CC ENAMEL ENL
1	30 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA shade indicator stick A3

VITAVM®CC VITA SYSTEM 3D-MASTER®



VITAVM®CC BASIC KIT 16-COLOR VITA SYSTEM 3D-MASTER®

Basic kit for BASIC layering

CVCCBK16

Quantity	Content	Material
16	30 g	VITA VM CC BASE DENTINE 0M1 - 5M1
2	30 g	VITA VM CC ENAMEL ENL, END
1	30 g	VITA VM CC WINDOW WIN
1	100 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITAVM®CC BASIC KIT 5-COLOR VITA SYSTEM 3D-MASTER®

Small basic kit for BASIC layering

CVCCBK5

Quantity	Content	Material
5	30 g	VITA VM CC BASE DENTINE 1M1, 2M2, 3M2, 3R2.5, 4M2
2	30 g	VITA VM CC ENAMEL ENL, END
1	30 g	VITA VM CC WINDOW WIN
1	30 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA Toothguide 3D-Master with Bleached Shade Guide



VITA VM CC ONE COLOR KIT 2M2

One-shade kit for initial tests

CVCCOCK2M2

Quantity	Content	Material
1	30 g	VITA VM CC BASE DENTINE 2M2
1	30 g	VITA VM CC ENAMEL ENL
1	30 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA shade indicator stick 2M2

VITA VM CC ONE COLOR KIT 3M2

One-shade kit for initial tests

CVCCOCK3M2

Quantity	Content	Material
1	30 g	VITA VM CC BASE DENTINE 3M2
1	30 g	VITA VM CC ENAMEL ENL
1	30 ml	VITA VM CC LIQUID
1		Measuring spoon, acrylic
1		VITA shade indicator stick 3M2

VITAVM®CC additional powder

compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades



VITAVM®CC PROFESSIONAL KIT

Additional materials for achieving special shade effects and individualizing the restoration

CVCCPK

Quantity	Content	Material
5	30 g	VITA VM CC EFFECT LINER EL2 - EL6

VITAVM® CC products

VITAVM® CC VITA classical A1-D4®

VITA VM CC BASE DENTINE, 30 g / 100 g

Description	Powder no.	Prod. no. 30 g	Prod. no. 100 g
A1	ENL	C4743130	C47431100
A2	ENL	C4743230	C47432100
A3	ENL	C4743330	C47433100
A3.5	END	C4743430	C47434100
A4	END	C4743530	C47435100
B2	END	C4743730	C47437100
B3	END	C4743830	C47438100
B4	END	C4743930	-
C1	END	C4744030	-
C2	END	C4744130	C47441100
C3	ENL	C4744230	C47442100
C4	ENL	C4744330	-
D2	END	C4744430	-
D3	END	C4744530	C47445100
D4	END	C4744630	-

VITAVM® CC VITA SYSTEM 3D-MASTER®

VITA VM CC BASE DENTINE, 30 g / 100 g

Description	Powder no.	Prod. no. 30 g	Prod. no. 100 g
0M1	ENL	C4703130	C47031100
1M1	ENL	C4703430	C47034100
2L1.5	ENL	C4703630	C47036100
2M1	ENL	C4703830	C47038100
2M2	ENL	C4703930	C47039100
2M3	ENL	C4704030	C47040100
3L1.5	ENL	C4704330	C47043100
3L2.5	ENL	C4704430	C47044100
3M1	ENL	C4704530	C47045100
3M2	ENL	C4704630	C47046100
3M3	ENL	C4704730	C47047100
3R1.5	ENL	C4704830	C47048100
3R2.5	ENL	C4704930	C47049100
4L1.5	END	C4705030	C47050100
4M2	END	C4705330	C47053100
5M1	END	C4705730	C47057100

VITAVM® CC additional powder

compatible with VITA SYSTEM 3D-MASTER and VITA classical A1-D4 shades

VITA VM CC EFFECT LINER, 30 g

Description	Prod. no. 30 g
EL2	C4714230
EL3	C4714330
EL4	C4714430
EL5	C4714530
EL6	C4714630

VITA VM CC WINDOW, 30 g

Description	Prod. no. 30 g
WIN	C4718130

VITA VM CC ENAMEL, 30 g / 100 g

Description	Prod. no. 30 g	Prod. no. 100 g
ENL	C4719130	C47191100
END	C4719230	C47192100

Description	Prod. no.
VITA VM CC LIQUID, 100 ml	CVCCCL100
VITA VM CC LIQUID, 30 ml	CVCCCL30

Accessories for Veneering Materials



Product description

The accessories for veneering materials make it easier to process these materials.

Accessories for Veneering Materials products

Accessories for Veneering Materials

Description		Prod. no.
Firing tray G (pack of 2 trays and 20 pins) 56 x 45 x 6 mm, with 52 holes for pins		B009N
Fibrous pad firing support, 2 pcs.		B017N
Pins for firing tray G (pack of 50)		B023
VITA glass spatula for stains		B055
VITA Karat diamond polishing set Contains 5 g diamond polishing paste, 20 diamond felt wheels, Ø 12 mm and 1 nickel-plated mandrel		B068
VITA Karat diamond polishing paste, 5 g		B069
VITA Karat diamond felt wheels, 20 pcs.		B070
VITA Modisol® Insulating pen with 2 tips (ceramic/plaster)		B075N
Brush No. 3 layering brush, red marten		B087
Brush No. E 0 staining		B093
Firing tray W for crowns and bridges		B201N
Platinum pins for firing tray W, 6 pcs.		B207N
Mixing spatula, Acrylic for PASTE OPAQUE		B227
VITA Silver Test Set To test the temperature inside ceramic furnaces, contains 6 ceramic sockets and 3 silver rods each 70 mm long and 1.5 mm Ø		B230
VITA INTERNO mixing palette, acrylic (Ø 144 mm, with indentations) Deckel ohne Bedruckung		B252
VITA AKZENT Plus mixing plate, crystal glass (70 x 95 mm)		B255V1
Shade indicator for VITA TITANKERAMIK 3D-MASTER Opaque, Opaque Dentine, Dentine, Enamel, Translucent, Window		B264
Shade guide for VITA TITANKERAMIK INTENSIVE / GINGIVA IN1-IN10, G1-G2		B274IN

Description		Prod. no.
Shade guide for TITANKERAMIK MARGIN / COLOR OPAQUE MARN, MAR1-MAR6, CO1-CO8		B274M
Shade guide for TITANKERAMIK TRANSLUCENT / CERVICAL T1-T6, WIN, CE1-CE4		B274T
Shade indicator for VITA TITANKERAMIK classical Opaque, Opaque Dentine, Dentine, Enamel, Translucent, Window		B288
Flat brush, for PASTE OPAQUE		B297
Stain brush tips for a & e brush no. 3/0, 5 pcs.		B303
Brush No. 8 a & e		B304
Tips for a & e brush no. 8, 2 pcs.		B305
Shade indicator for VITA VM 9 3D-MASTER Base Dentine, Transpa Dentine, Enamel, Chroma Plus		B330
Shade indicator for VITA VM 9 classical Base Dentine, Transpa Dentine, Enamel, Chroma Plus		B331
Shade indicator for VITA VM 13 3D-MASTER Opaque, Sun Dentine, Base Dentine, Transpa Dentine, Enamel, Chroma Plus, Neutral, Window, Corrective, Color Opaque		B340
VITA Shade indicator for VITA VM 13 classical Opaque, Base Dentine, Transpa Dentine, Enamel, Chroma Plus, Neutral, Window, Corrective, Color Opaque		B341
Shade indicator for VITA VMK MASTER classical Opaque, Opaque Dentine, Dentine, Enamel, Wash Opaque, Translucent T4, WIN		B342
Shade indicator for VITA VMK MASTER 3D-MASTER Opaque, Opaque Dentine, Dentine, Enamel, Wash Opaque, Translucent T4, WIN		B343
Shade guide for VITA VM EFFECT CHROMA EC1-EC11		B414
Shade guide for VITA VM EFFECT ENAMEL / EFFECT PEARL EE1-EE11, EP1-EP3		B416
Shade guide for VITA VM PROFESSIONAL KIT SMALL EC1, EC4, EC6, EC8, EC9, MM2, EP1, EO2, EE1, EE3, EE7, EE8, EE9, EE10, EE11		B419
Shade guide for VITA VM INDIVIDUAL For creating individual shade indicator plates		B423
Shade guide for VITA VMK MASTER TRANSLUCENT T1-T8, OT1, PLT1		B431
Shade guide for VITA VMK MASTER MARGIN / GINGIVA MN, M1-M5, G1-G5, GOL, GOD		B432
Shade guide for VITA VMK MASTER DENTINE MODIFIER / LUMINARY DM1-DM8, LM1-LM6		B433
Shade guide for VITA VMK MASTER STANDARD N1-N5, CE1-CE3, COR1, COR2		B434
Shade guide for VITA VM 9 ADD-ON ADD1 - ADD8		B435

Description		Prod. no.
VITA AKZENT Plus Shade sample disc		B436
VITA AKZENT Plus mixing plate, acrylic (Ø 144 mm, with indentations)		B437
Shade guide for VITA VM GINGIVA / MARGIN G1-G5, GOL, GOD, M1-M9		B439
Shade guide for VITA VM EFFECT LINER / EFFECT OPAL / MAMELON EL1-EL6, E01-E05, MM1-MM3		B440
Shade guide for VITA VM 11 No. 1 MM1, MM3, SD1, SD2, SD3, EC1, EC5, EC11		B442
Shade guide for VITA VM 11 No. 2 E01, E02, E03, E05, EP1, EE1, EE3, EE5, EE7, EE8, EE9, EE11		B443
VITA INTERNO Shade sample disc		B444

Description	Prod. no.
VITA AKZENT Plus CHROMA STAINS, color sample classical	B455
VITA AKZENT Plus CHROMA STAINS, color sample 3D-MASTER	B458

Liquids for Veneering Materials



Product description

Modeling liquids for outstanding stability when layering dental ceramics.

Liquids for Veneering Materials products

Liquids for Veneering Materials

Description		Prod. no.
VITA MODELLING LIQUID 30M, 250 ml For mixing dentine and enamel materials when modeling VITA porcelains for long periods with firing temperatures of over 900°C		B30M250
VITA MODELLING LIQUID 30M, 50 ml For mixing dentine and enamel materials when modeling VITA porcelains for long periods with firing temperatures of over 900°C		B30M50
VITA HIGH SILVER MODELLING LIQUID, 50 ml Anti-greening liquid for high silver content alloys with a silver content of more than 30 %		BHSM50
VITA LUMEX AC MODELLING LIQUID, 250 ml		BLML250
VITA LUMEX AC MODELLING LIQUID, 50 ml		BLML50
VITA MODELLING FLUID, 250 ml Standard modeling liquid For mixing VITA VM 13, VMK Master- and TITANKERAMIK powders		BMF250
VITA MODELLING FLUID, 50 ml Standard modeling liquid For mixing VITA VM 13, VMK Master- and TITANKERAMIK powders		BMF50
VITA MODELLING FLUID RS, 250 ml Special modelling liquid for modelling over an extended period For mixing any dentin, incisal or additional materials. Its smooth consistency allows extended and wet processing while also ensuring good stability. This fluid is perfectly suited for use in larger restorations or multi-unit bridges.		BMFRS250
VITA MODELLING FLUID RS, 50 ml Special modelling liquid for modelling over an extended period For mixing any dentin, incisal or additional materials. Its smooth consistency allows extended and wet processing while also ensuring good stability. This fluid is perfectly suited for use in larger restorations or multi-unit bridges.		BMFRS50
VITA VM MODELLING LIQUID, 250 ml Standard modeling liquid For mixing VITA VM 7, VM 9, VM 13 BASE DENTINE, TRANSPA DENTINE, ENAMEL and EFFECT materials		BML250N
VITA VM MODELLING LIQUID, 50 ml Standard modeling liquid For mixing VITA VM 7, VM 9, VM 13 BASE DENTINE, TRANSPA DENTINE, ENAMEL and EFFECT materials		BML50N
VITA OPAQUE FLUID, 250 ml For mixing opaque and bonder For mixing VITA VMK Master- and TITANKERAMIK powders		BOF250

Description		Prod. no.
VITA OPAQUE FLUID, 50 ml For mixing opaque and bonder For mixing VITA VMK Master- and TITANKERAMIK powders		BOF50
VITA VM PASTE FLUID, 15 ml To change the consistency of OPAQUE PASTE and COLOR OPAQUE PASTE For mixing VITA VM 13 and VM 15 OPAQUE pastes		BPF15V1
VITA PASTE OPAQUE LIQUID, 15 ml To change the consistency of VITA paste opaques and paste bonders		BPOL15
VITA die separating liquid, 16 ml Separator (ceramic/plaster) for use with VITA ceramics (except TITANKERAMIK)		BSI16
VITA LOW FUSING MODELLING LIQUID, 50 ml Special mixing fluid for low-melting ceramic Liquid for mixing the VITA VM 11 TRANSPA DENTINE, SUN DENTINE, ENAMEL and additional materials		BVMLF50
VITA VM OPAQUE FLUID, 250 ml For mixing opaque powder For mixing VITA VM 13 and VM 15 OPAQUE materials		BVMOF250
VITA VM OPAQUE FLUID, 50 ml For mixing opaque powder For mixing VITA VM 13 and VM 15 OPAQUE materials		BVMOF50

Liquids for Veneering resin

Description		Prod. no.
VITA AKZENT LC CLEANER, 25 ml		CALCC25
VITA VM LC CLEANER, 50 ml Cleaning solution		CVMLC50
VITA VM LC GEL, 20 ml Inhibition layer blocker		CVMLG20
VITA VM LC MODELLING LIQUID, 10 ml Modeling liquid		CVMLM10
Modeling liquid Modeling liquid		CVMLM30
VITA VM LC SEPARATOR, 30 ml Modeling separator		CVMLS30

Accessories for Veneering resin



Product description

The accessories for veneering acrylics make it easier to process these materials.

Accessories for Veneering resin products

Accessories for Veneering resin



Description		Prod. no.
Measuring spoon, acrylic		C006
VITA Mixing palette black porcelain with lid, 8,5 x 11 cm		C014
Measuring spoon, 0,1 ml / Mixing spatula 2-sided		C017
VITA Mixing palette white porcelain with lid, 8,5 x 11 cm		C018
White brush-holder		FCPHW
Disposable brush heads, 50 pcs.		FCPS

VITA SPRAY-ON



Product description

VITA SPRAY-ON system is for spraying dental ceramic opaque materials onto metal frameworks.

A compressed air spray gun is used to spray on a uniformly thin covering layer of opaque material. Materials which have already been mixed can be kept in the glass bottles supplied with the system until they have been completely used up.

Adaptors of varying diameters are supplied for connection to suction units.

Features/Benefits

- Extremely time-saving
- Produces a uniform and extremely thin covering
- Achieves optimal bonding with metal framework
- Considerably reduces costs

VITA SPRAY-ON kits

VITA SPRAY-ON

VITA SPAY-ON Sortiment / Assortment B101V3

Quantity	Content	Material
13		VITA SPRAY-ON Glass bottle empty
1		VITA SPRAY-ON Glass bottle empty, OP0
1		VITA SPRAY-ON Glass bottle empty, OP1
1		VITA SPRAY-ON Glass bottle empty, OP2
1		VITA SPRAY-ON Glass bottle empty, OP3
1		VITA SPRAY-ON Glass bottle empty, OP4
1		VITA SPRAY-ON Glass bottle empty, OP5
1		VITA SPRAY-ON Glass bottle empty, WO
1		Sticker VITA classical A1-D4 (16 pieces)
1	6 pcs.	VITA SPRAY-ON Adapter lid with sealing washers
1		VITA SPRAY-ON Spray gun with protective cap
1		VITA SPRAY-ON Replacement needle for spray gun
1	2 pcs.	VITA SPRAY-ON Teflon needle bearings
1		VITA glass spatula for stains
1		VITA SPRAY-ON Holder for two spray guns
1		VITA SPRAY-ON Air hose with woven fabric inlay
1		VITA SPRAY-ON Spiral connector
1	1000 ml	VITA SPRAY-ON Liquid
1	250 ml	VITA SPRAY-ON Indikator-Liquid
1	250 ml	VITA SPRAY-ON Dispensing bottle empty
1		VITA SPRAY-ON Permanent marker
1		VITA SPRAY-ON Pressure regulator with wall mounting

VITA SPRAY-ON products**VITA SPRAY-ON**

Description	Prod. no.
VITA SPRAY-ON Glass bottle empty	B103
VITA SPRAY-ON Adapter lid with sealing washers, 6 pcs.	B104
VITA SPRAY-ON Spray gun with protective cap	B105N
VITA SPRAY-ON Replacement needle for spray gun	B106
VITA SPRAY-ON Pressure regulator with wall mounting	B107N
VITA SPRAY-ON Holder for two spray guns	B108
VITA SPRAY-ON Teflon needle bearings, 2 pcs.	B109N
VITA SPRAY-ON Air hose with woven fabric inlay	B110
VITA SPRAY-ON Spiral connector	B111
VITA SPRAY-ON Liquid, 1000 ml	B112
VITA SPRAY-ON Indikator-Liquid, 250 ml	B113
VITA SPRAY-ON Dispensing bottle empty, 250 ml	B114
VITA SPRAY-ON Protective cap for spray gun	B115N
VITA SPRAY-ON Permanent marker	B116
VITA SPRAY-ON Sealing washers for adapter cable, 10 pcs.	B210

VITA Modulbox

Plastic cabinet storage system



Product description


The VITA MODULBOX is a modular A3-sized system for long-term, clearly-arranged and material-appropriate storage of VITA products in laboratories and dental practices.

See also the vario tooth cabinet.

VITA Modulbox products

VITA Modulbox

Description		Prod. no.
VITA Modulbox 20 mm drawer with removable section and 2 labelling fields		ZMB20
VITA Modulbox 55 mm drawer with removable section and 2 labelling fields		ZMB55
VITA Modulbox 80 mm drawer with removable section and 2 labelling fields		ZMB80
VITA Modulbox top and bottom sections without removable section, incl. Accessories		ZMBBD
VITABLOCS for CEREC Storage Box with 1 drawer 80 mm with insert No. 7		ZMBCER1
VITA Modulbox Storage Box No. 1 with 1 drawer 55 mm		ZMBK1
VITA Modulbox Storage Box No. 2 with 2 drawers 55 mm		ZMBK2
VITA Modulbox Storage Box No. 3 with 1 drawer 80 mm		ZMBK3
VITA Modulbox Storage Box No. 4 with 2 drawers 80 mm		ZMBK4
VITA Modulbox Storage Box No. 5 with 2 drawers, each 1 x 55 mm and 80 mm		ZMBK5
Tooth cabinet S10 with 10 drawers: External dimensions: 482 x 345 x 350 mm		ZMBS10
Tooth cabinet S12 with 12 drawers: External dimensions: 482 x 345 x 420 mm		ZMBS12
Tooth cabinet S2 with 2 drawers: External dimensions: 482 x 345 x 70 mm		ZMBS2
Tooth cabinet S4 with 4 drawers: External dimensions: 482 x 345 x 140 mm		ZMBS4
Tooth cabinet S6 with 6 drawers: External dimensions: 482 x 345 x 210 mm		ZMBS6
Tooth cabinet S8 with 8 drawers: External dimensions: 482 x 345 x 280 mm		ZMBS8
VITA Modulbox intermediate sections without removable section, incl. 3 mountings		ZMBZ

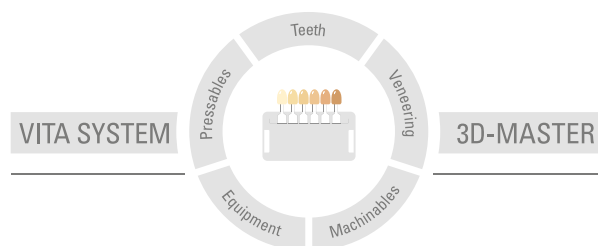

This information issued: 09-2021

Our products were developed specifically for dental use and must be processed according to the instructions for use. The information/illustrations contained herein do not constitute an assurance of properties and are not binding.

VITA

VITA Zahnfabrik H. Rauter GmbH & Co.KG
Postfach 1338 · D-79704 Bad Säckingen · Germany
Tel. +49 (0) 7761/562-0 · Fax +49 (0) 7761/562-299
Hotline: Tel. +49 (0) 7761/562-222 · Fax +49 (0) 7761/562-446
www.vita-zahnfabrik.com · info@vita-zahnfabrik.com
 facebook.com/vita.zahnfabrik