
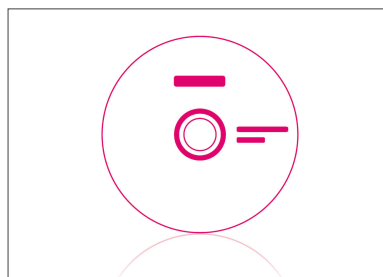
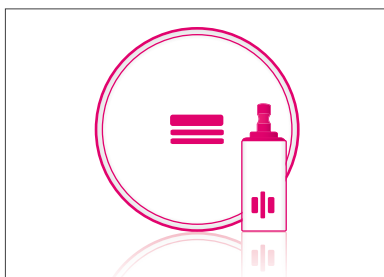
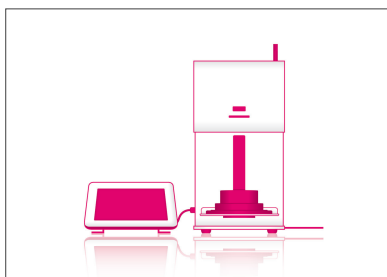
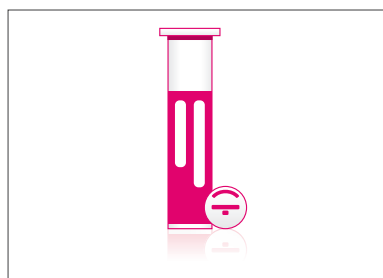


# VITA Compendium September 2021

## Luting

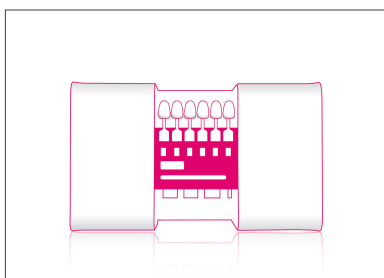

VITA – perfect match.

**VITA**



VITA – perfect match.

**VITA**

## VITA ADIVA® FULL-ADHESIVE

The full-adhesive luting system with perfectly matched components for conditioning tooth substance and restorative material



### Product description

The full-adhesive luting system with perfectly matched components for conditioning tooth substance and restorative material.

Suitable for luting of all indirect restorations made of feldspar, glass and hybrid ceramics, zirconia, composite and metal.

### Features/Benefits

- Luting numerous restorative materials with material-specific, matched coupling agents in a dependable way.
- Process-proof luting thanks to a clearly structured luting set with systematically arranged and color-coded components.
- All components are easily accessible thanks to the handy two-piece set in a standard tray.

### Content

#### VITA ADIVA® TOOTH-ETCH

is a 37% phosphoric acid gel for etching tooth enamel for adhesive bonding of restorations to natural tooth substance. Also suitable for the total-etch and selective-etch techniques.

#### VITA ADIVA® T-BOND I + II

is a dual-curing bonder system. VITA ADIVA T-BOND provides an excellent bond between the luting composite and etched enamel, as well as dentin.

#### VITA ADIVA® CERA-ETCH

is a 5% hydrofluoric acid gel used for etching restorations made of feldspar, glass and hybrid ceramics, to achieve a micromechanical bond between the restoration and luting composite. Intended for extraoral use only!

#### VITA ADIVA® C-PRIME

is a coupling agent that creates a chemical bond between restorations made of feldspar, glass and hybrid ceramics and luting composite.

#### VITA ADIVA® ZR-PRIME

is a coupling agent that creates a chemical bond between sintered zirconia restorations and luting composite.

#### VITA ADIVA® M-PRIME

is a coupling agent that creates a chemical bond between sandblasted restorations made of titanium, non-precious metal alloys, precious metal alloys and luting composite.

#### VITA ADIVA® F-CEM

is a dual-curing, full-adhesive composite that's suitable for permanent luting of restorations made of feldspar, glass and hybrid ceramics, zirconia, composite and metal.

#### VITA ADIVA® OXY-PREVENT

is a neutral-colored glycerine gel used during luting to prevent the formation of an oxygen inhibition layer. Can also be used as a try-in paste.

## VITA ADIVA® FULL-ADHESIVE kits

### VITA ADIVA® FULL-ADHESIVE LUTING SET

#### VITA ADIVA® FULL-ADHESIVE LUTING SET FAFAS

Quantity	Content	Material
1	5 ml	VITA ADIVA F-CEM A2 Universal
1	3 ml	VITA ADIVA OXY-PREVENT
1	10 pcs.	Disposable application cannulae, green for OXY-PREVENT
1	3 ml	VITA ADIVA TOOTH-ETCH
1	for TOOTH	Disposable application needles blue
1	3 ml	VITA ADIVA CERA-ETCH
1	20 pcs	Disposable application cannulae red for CERA ETCH
1	5 ml	VITA ADIVA T-BOND I
1	5 ml	VITA ADIVA T-BOND II
1	5 ml	VITA ADIVA C-PRIME
1	5 ml	VITA ADIVA ZR-PRIME
1	5 ml	VITA ADIVA M-PRIME
1	25 pcs.	VITA ADIVA Microbrush
1	12 pcs.	VITA ADIVA Mixing-Trays
1	15 pcs.	VITA ADIVA T-Mixer 1:1

## VITA ADIVA® FULL-ADHESIVE LUTING VERSIONS

VITA ADIVA F-CEM, 5 ml

Description	Prod. no. 5 ml
A2 Universal	FAFCEMA25
A3	FAFCEMA35
Translucent	FAFCEMT5
White opaque	FAFCEMW05

## VITA ADIVA® COMPONENTS

Description	Prod. no.
VITA ADIVA CERA-ETCH, 3 ml 5% hydrofluoric acid gel	FACE3
VITA ADIVA CERA-ETCH, 6 ml 5% hydrofluoric acid gel	FACE6
VITA ADIVA C-PRIME, 5 ml Ceramic primer	FACP5
VITA ADIVA M-PRIME, 5 ml Metal / Alloy primer	FAMP5
VITA ADIVA OXY-PREVENT, 3 ml Oxygen protection gel	FAOP3
VITA ADIVA T-BOND SET (only for VITA ADIVA FULL-ADHESIVE) Dual cure adhesive system	FATBS
VITA ADIVA TOOTH-ETCH, 3 ml (only for VITA ADIVA FULL-ADHESIVE) 37% phosphoric acid gel	FATE3
VITA ADIVA ZR-PRIME, 5 ml Zirconia primer	FAZRP5

**VITA ADIVA® ACCESSORIES**

Description	Prod. no.
VITA ADIVA Mixing-Trays, 25 pcs.	FAAT25
VITA ADIVA Microbrush, 25 pcs.	FAMB25
VITA ADIVA T-Mixer 1:1, 15 pcs. (only for VITA ADIVA FULL-ADHESIVE)	FATM11
VITA ADIVA T-Mixer 4:1, 15 pcs. (only for VITA ADIVA SELF-ADHESIVE)	FATM41

## VITA ADIVA® SELF-ADHESIVE

The self-adhesive luting system with perfectly matched components for conditioning the restorative material



### Product description

The self-adhesive luting system with perfectly matched components for conditioning the restorative material.

Suitable for luting of all indirect restorations made of feldspar, glass and hybrid ceramics, zirconia, composite and metal.

### Features/Benefits

- Luting numerous restorative materials with material-specific, matched coupling agents in a dependable way.
- Process-proof luting thanks to a clearly structured luting set with systematically arranged and color-coded components.
- All components are easily accessible thanks to the handy two-piece set in a standard tray.

### Content

#### VITA ADIVA® CERA-ETCH

is a 5 % hydrofluoric acid gel used for etching restorations made of feldspar, glass and hybrid ceramics, to achieve a micromechanical bond between the restoration and luting composite. Intended for extraoral use only!

#### VITA ADIVA® C-PRIME

is a coupling agent that creates a chemical bond between restorations made of feldspar, glass and hybrid ceramics and luting composite.

#### VITA ADIVA® ZR-PRIME

is a coupling agent that creates a chemical bond between sintered zirconia restorations and luting composite.

#### VITA ADIVA® M-PRIME

is a coupling agent that creates a chemical bond between sandblasted restorations made of titanium, non-precious metal alloys, precious metal alloys and luting composite.

#### VITA ADIVA® S-CEM

is a dual-curing, self-adhesive composite that's suitable for permanent luting of restorations made of feldspar, glass and hybrid ceramics, zirconia, composite and metal.

#### VITA ADIVA® OXY-PREVENT

is a neutral-colored glycerine gel used during luting to prevent the formation of an oxygen inhibition layer. Can also be used as a try-in paste.



## VITA ADIVA® SELF-ADHESIVE kits

## VITA ADIVA® SELF-ADHESIVE LUTING SET

### VITA ADIVA® SELF-ADHESIVE LUTING SET FASAS

Quantity	Content	Material
1	5 ml	VITA ADIVA S-CEM A2 Universal
1	3 ml	VITA ADIVA OXY-PREVENT
1	10 pcs.	Disposable application cannulae, green for OXY-PREVENT
1	3 ml	VITA ADIVA CERA-ETCH
1	20 pcs	Disposable application cannulae red for CERA ETCH
1	5 ml	VITA ADIVA C-PRIME
1	5 ml	VITA ADIVA ZR-PRIME
1	5 ml	VITA ADIVA M-PRIME
1	25 pcs.	VITA ADIVA Microbrush
1	12 pcs.	VITA ADIVA Mixing-Trays
1	15 pcs.	VITA ADIVA T-Mixer 4:1

## VITA ADIVA® SELF-ADHESIVE LUTING VERSIONS

VITA ADIVA S-CEM, 5 ml

Description	Prod. no. 5 ml
A2 Universal	FASCEMA25
Translucent	FASCEMT5

## VITA ADIVA® COMPONENTS

Description	Prod. no.
VITA ADIVA CERA-ETCH, 3 ml 5% hydrofluoric acid gel	FACE3
VITA ADIVA CERA-ETCH, 6 ml 5% hydrofluoric acid gel	FACE6
VITA ADIVA C-PRIME, 5 ml Ceramic primer	FACP5
VITA ADIVA M-PRIME, 5 ml Metal / Alloy primer	FAMP5
VITA ADIVA OXY-PREVENT, 3 ml Oxygen protection gel	FAOP3
VITA ADIVA T-BOND SET (only for VITA ADIVA FULL-ADHESIVE) Dual cure adhesive system	FATBS
VITA ADIVA TOOTH-ETCH, 3 ml (only for VITA ADIVA FULL-ADHESIVE) 37% phosphoric acid gel	FATE3
VITA ADIVA ZR-PRIME, 5 ml Zirconia primer	FAZRP5

## VITA ADIVA® ACCESSORIES

Description	Prod. no.
VITA ADIVA Mixing-Trays, 25 pcs.	FAAT25
VITA ADIVA Microbrush, 25 pcs.	FAMB25
VITA ADIVA T-Mixer 1:1, 15 pcs. (only for VITA ADIVA FULL-ADHESIVE)	FATM11
VITA ADIVA T-Mixer 4:1, 15 pcs. (only for VITA ADIVA SELF-ADHESIVE)	FATM41

## VITA ADIVA® TE-CEM



### Product description

VITA ADIVA TE-CEM is a dual-hardening composite cement for temporary bonding.

It is available in the Translucent shade.

## VITA ADIVA® TE-CEM products

### VITA ADIVA® TE-CEM

VITA ADIVA TE-CEM, 5 ml

Description	Prod. no. 5 ml
Translucent	FATECEM5

## VITA ADIVA® IA-CEM

Ultra opaque luting solution for implant restorations



### Product description

VITA ADIVA IA-CEM is an ultra opaque, strongly-masking, dual-curing luting composite for implant prosthetics. It is used for permanent adhesive bonding of (abutment) crowns and mesostructures on adhesive titanium bases and individual titanium and zirconia abutments. There are material-specific, matched bonding agents for the luting composite to ensure durable bonding of the restoration and the abutment.

## VITA ADIVA® IA-CEM kits

### VITA ADIVA® IA-CEM

#### VITA ADIVA® IA-CEM SET

#### VITA ADIVA® IA-CEM SET

#### FAIACEMSET

Quantity	Content	Material
1	5 ml	VITA ADIVA IA-CEM Ultra opaque
1	15 pcs.	VITA ADIVA T-Mixer 1:1
1	6 ml	VITA ADIVA CERA-ETCH
1	5 ml	VITA ADIVA C-PRIME
1	5 ml	VITA ADIVA ZR-PRIME
1	25 pcs.	VITA ADIVA Microbrush
1	12 pcs.	VITA ADIVA Mixing-Trays




**VITA ADIVA® IA-CEM products**

**VITA ADIVA® IA-CEM**

VITA ADIVA IA-CEM, 5 ml

Description	Prod. no. 5 ml
Ultra opaque	FAIACEMU05

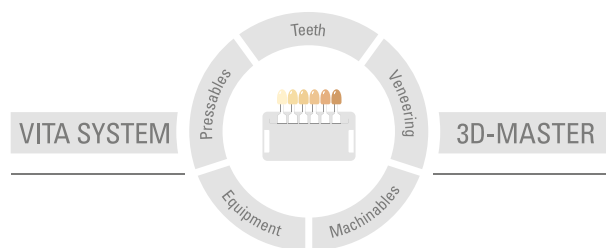

This information issued: 09-2021

Our products were developed specifically for dental use and must be processed according to the instructions for use. The information/illustrations contained herein do not constitute an assurance of properties and are not binding.

# VITA

VITA Zahnfabrik H. Rauter GmbH & Co.KG  
Postfach 1338 · D-79704 Bad Säckingen · Germany  
Tel. +49 (0) 7761/562-0 · Fax +49 (0) 7761/562-299  
Hotline: Tel. +49 (0) 7761/562-222 · Fax +49 (0) 7761/562-446  
[www.vita-zahnfabrik.com](http://www.vita-zahnfabrik.com) · [info@vita-zahnfabrik.com](mailto:info@vita-zahnfabrik.com)  
 [facebook.com/vita.zahnfabrik](https://facebook.com/vita.zahnfabrik)